# GEMINI **EVO**<sup>TM</sup>

810 + 980 DIODE LASER

---

USER MANUAL

## WARNINGS & CAUTIONS

Failure to comply with precautions and warnings described in this User Manual may lead to exposure to dangerous optical radiation sources. Please comply with all safety instructions and warnings.

### CAUTION:

- Read these instructions carefully prior to using your Gemini EVO 810+980 soft tissue diode laser.
- Ensure that all users are properly trained prior to use. Consult your distributor for training recommendations. Mandatory training on the Gemini EVO laser is done via this manual.
- Do not modify this equipment without authorization of the manufacturer.
- Laser fume and/or plume may contain viable tissue particles.
- Always wind the fiber optic cable in a clockwise manner around the fiber wrap to avoid fiber breakage.
- Do not use in the presence of combustible or combustion supporting gases.
- Always test activate the device outside the mouth before using on a patient.
- This unit has been designed and tested to meet the requirements of electromagnetic, electrostatic, and radio frequency interference standards. However, the possibility of electromagnetic or other interference may still exist. Relocating the device may help to eliminate the interference.
- Medical electrical equipment needs special precautions regarding electromagnetic compatibility (EMC) and needs to be installed and put into service according to the EMC information provided in this manual.
- Periodically inspect laser eyewear for pitting and cracking.
- In the event abnormal performance is observed, discontinue use and follow instructions in the service and troubleshooting section of this manual.

Safety is paramount when using any energy-based surgical instrument and your office should implement a safety program for the Gemini EVO 810+980 soft tissue diode laser. If your office does not already have a safety officer, one should be appointed to be responsible for understanding proper use, safe operation, and maintenance of the Gemini EVO laser system. Their duties should include training office personnel in all aspects of system safety and management of the Gemini EVO laser and all accessories. / **ADDITIONAL TROUBLESHOOTING:** For additional troubleshooting questions and tips call 1.801.553.4574. To check for the latest software updates, download the Gemini EVO App in the iOS or Android web store.

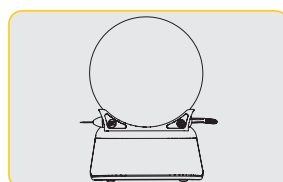
### WARNING:

- Visible and Invisible Laser Radiation – Avoid eye or skin exposure to direct or scattered radiation. Class IV Laser Product
- Laser Safety Eye Protection **MUST BE WORN** by the Operator, Patient, Assistant, and anyone present when the laser is activated. Eye Protection must conform to Specification DIN EN207 Annex II of the Directive 89/686/EEC with wavelength protection of 810nm- 980nm, and  $\pm 10$ nm of OD 5+.
- Never direct or point the beam at a person's eyes.
- Do not look directly into the beam or at specular reflections.
- Do not aim the laser at metallic or reflective surfaces, such as surgical instruments or dental mirrors. If aimed directly at these surfaces the laser beam will reflect and create a potential hazard.
- Never operate the laser without a fiber tip attached.
- Laser aperture at the end of the handpiece.
- Laser aperture warning label affixed to system handpiece.
- Always place the system into STBY mode when leaving the Gemini EVO 810+980 soft tissue laser unattended for a few minutes or between patients.
- Use of controls or adjustments or performance of procedures other than those specified herein may result in hazardous radiation exposure.
- Do not open unit housing at any time. Danger from optical radiation may exist.
- The use of accessories, other than those specified, except those supplied or sold by Ultradent Products, Inc., as replacement parts for internal or external components, may result in increased EMISSIONS or decreased IMMUNITY of the Gemini EVO 810+980 soft tissue laser.

## WHAT IS IN THE BOX


The Gemini EVO 810+980 soft tissue laser includes the following.


Laser Unit

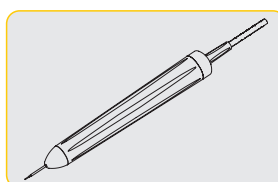
Activation Pedal with 2 AA Batteries

Fiber Delivery System

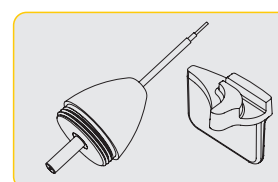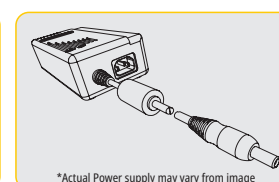Pre-Initiated Disposable Tips (10)  
Bending tool included in the tip box


\*Actual Power supply may vary from image  
\*DC Power Supply with Plug Adapters

\*Actual eyewear may vary from image

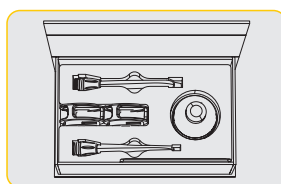\*Protective Eyewear (3)

3 mm, 7 mm, and 25 mm PBM  
adapters

- Laser Warning Sign
- USB Cable
- User Manual with Warranty Information

NOTE: The laser ships with the lithium ion battery and fiber delivery system already installed

NOTE: Use proper care when transporting the unit

ALSO INCLUDED: Laser Warning Sign and User's Manual

**WARNING: No modification of this equipment is allowed**

## UNPACKING INSTRUCTIONS

A manufacturer or dealer representative can provide assistance when you are ready to remove the laser from its shipping container. Please do not attempt to unpack the Gemini EVO 810+980 soft tissue laser or install the system without reading this manual first. If you are unsure about any aspect of the assembly, call your customer service representative or dealer for assistance.

## SHIPPING CONTAINER INFORMATION

The shipping container you received with your Gemini EVO 810+980 soft tissue laser was specially designed to safely transport the device. In the unlikely event that you need to return the laser for service or repair, please retain the original shipping container.

# TABLE OF CONTENTS

---

## OVERVIEW

Display .....	3
Guided Touch Interface "GTI" .....	3
Laser Unit .....	3
Activation Pedal .....	3
Fiber Delivery System .....	3
Quick Start .....	4

---

## CONTROLS - OPERATION & USE

01 - Electronic Key Passcode .....	4
02 - Selecting a Wavelength .....	4
03 - Activation Pedal Connection .....	4
04 - Manual Power Adjustment .....	4
05 - Laser STBY and Active Modes .....	4-5
06 - Preset Procedures Settings and Customization .....	5
07 - Tip Illumination .....	5
08 - Haptic Sense "HS" - Activation Pedal and PBM adjustments .....	5-6
09 - Sound .....	6
10 - Aiming Light .....	6
11 - Photobiomodulation (PBM) .....	6-7
12 - Disposable Tip Operation .....	8
13 - Battery and Battery Level Indications .....	8
14 - Power Supply .....	8
15 - Fiber Wrapping .....	8-9
16 - Handpiece Magnet .....	9
17 - Operating Mode .....	9
18 - Emergency Stop .....	9
19 - Remote Interlock .....	9
20 - Adding New Activation Pedal and Pairing .....	9
21 - Transparent Electroluminescent Display .....	9
22 - Wi-Fi Connectivity and Dashboard .....	9-12
Cleaning and Sterilization Procedures .....	12
Procedural Recommendations .....	12-13
Facility & Environmental Considerations .....	13
General Safety Considerations .....	13-14
System Specifications .....	14-15
Service and Troubleshooting .....	15-16
Cybersecurity .....	17
Electromagnetic Environment Guidance .....	17-19
Warranty .....	19
Labeling .....	20

## OVERVIEW - DISPLAY

### SEE FIGURE 1.1

01 - ELECTRONIC PASSCODE	07 - POWER BAR	13 - BATTERY	19 - 980 ACTIVE/STBY
02 - PRESET PROCEDURES	08 - VOLUME INDICATION	14 - POWER SUPPLY	20 - VALUE TYPE INDICATOR
03 - PRESET CATEGORIES	09 - REMOTE INTERLOCK	15 - POWER INDICATOR	
04 - 810 ACTIVE/STBY	10 - SOFTWARE UPDATE	16 - SOFTWARE VERSION	
05 - 810 WAVELENGTH	11 - WI-FI STRENGTH	17 - AIMING LIGHT LEVEL	
06 - SOUND LEVEL	12 - BLUETOOTH	18 - 980 WAVELENGTH	

## OVERVIEW - GUIDED TOUCH INTERFACE "GTI"

### SEE FIGURE 1.2

01 - EMERGENCY STOP	08 - VOICE CONFIRMATION	15 - HAPTIC SENSE
02 - 810 WAVELENGTH	09 - TIP ILLUMINATION	16 - AVAILABLE UPDATES
03 - DUAL WAVELENGTH	10 - SOUND	17 - PRESET PROCEDURES
04 - 980 WAVELENGTH	11 - ACTIVE/STBY	18 - HOME
05 - ON/OFF	12 - UP/DOWN ARROWS	
06 - WI-FI INDICATOR	13 - LEFT/RIGHT ARROWS	
07 - MENU SELECTION	14 - AIMING LIGHT	

## OVERVIEW - LASER UNIT

### SEE FIGURE 1.3

01 - TRANSPARENT DISPLAY
02 - BASE
03 - STRAIN RELIEF
04 - ON/OFF
05 - GUIDED TOUCH INTERFACE
06 - EMERGENCY STOP

### SEE FIGURE 1.4

01 - MAGNETIC HANDPIECE HOLDER	05 - OPTICAL FIBER WITH PROTECTION TUBE
02 - LED TIP	06 - REMOTE INTERLOCK PORT
03 - OPTICAL FIBER	07 - EXTERNAL POWER SUPPLY PLUG
04 - SURGICAL HANDPIECE	08 - MINI USB (FACTORY USE ONLY)

## OVERVIEW - ACTIVATION PEDAL

### SEE FIGURE 1.5

01 - PEDAL STATUS INDICATORS	07 - BATTERY INDICATOR	09 - HAPTIC SENSE "HS" INDICATOR
02 - OVERSTEP / SHROUD PROTECTION	7.1 - 60 TO 100%	10 - BLUETOOTH (ICON TURNS BLUE WHEN CONNECTED)
03 - ACTIVATION AREA	7.2 - 25 TO 59%	
04 - BATTERY COMPARTMENT LID	7.3 - 0 TO 24%	
05 - AA BATTERIES (X2)	08 - LASER ACTIVE TURNS GREEN WHEN PEDAL IS PRESSED	
06 - LABEL / ANTI SKID AREAS		

## OVERVIEW - FIBER DELIVERY SYSTEM

### SEE FIGURE 1.6

The Fiber Delivery System is a unique and ergonomic optical cable that is non-detachable from the Laser Unit. The handpiece will require cleaning and sterilization after each patient treatment. Disposable tips are intended for single-use only and must be disposed of after each patient use.

01 - ANODIZED ALUMINUM HANDPIECE	06 - HAPTIC SENSE "HS" - ENABLED ONLY WITH PBM PROCEDURES
02 - ILLUMINATION TIP	07 - TURN COUNTER CLOCKWISE (REMOVE TO STERILIZE)
03 - OPTICAL FIBER	08 - FIBER STRAIN RELIEF
04 - BENDABLE STAINLESS STEEL TUBE	
05 - OPTICAL FERRULE (DO NOT TOUCH)	

# QUICK START

## 1. DOWNLOAD THE MOBILE APP (Figure: 2.1)

## 2. PLUG IN POWER SUPPLY

During initial setup use the AC/DC power supply for at least one hour to fully charge the battery. Plug the power supply into an AC outlet and connect to the corresponding connector on the rear of the system. (Figure: 2.2)

## 3. INSERT AA BATTERIES INTO ACTIVATION PEDAL

Install the (2) provided AA batteries into the wireless activation pedal. When replacing the AA batteries, we recommend an ALKALINE type battery. (Figure: 2.3)

## 4. TURN LASER UNIT ON

The universal ON/OFF button is a membrane switch that requires pressure in order to be activated. (Figure: 2.4)

## 5. ENTER ELECTRONIC KEY PASSCODE

Enter the electronic key passcode on the Guided Touch Interface using the Up and Down arrow keys. The security code sequence is UP, DOWN, UP, DOWN. A checkmark icon will appear when the correct key is input. (Figure: 2.5)

## 6. SELECT YOUR DESIRED WAVELENGTH

Select the desired laser wavelength on the Guided Touch Interface: 810 nm, 980 nm, or Dual Wavelength. (Figure: 2.6)

"Please select wavelength" ➡

## 7. SELECT YOUR DESIRED POWER SETTING

Select your desired power setting, then activate the laser. (Figure: 2.7)

# CONTROLS, OPERATION & USE

## 01 - ELECTRONIC KEY PASSCODE

The Gemini EVO 810+980 soft tissue laser is equipped with an electronic key passcode. When you turn the Laser Unit on, the passcode key screen will be displayed at the bottom center of the screen. The correct passcode sequence should be entered on the Guided Touch Interface: UP, DOWN, UP, DOWN. (Figure: 3.1, Options: 1, 2, 1, 2)

The Gemini EVO 810+980 soft tissue laser is equipped with Guided Touch Interface "GTI" which means only the icons that are relevant for a given procedure will be shown. When entering the electronic key passcode, only the UP and DOWN arrows will be shown as they are the only necessary icons to be touched while entering the passcode.

THE GUIDED TOUCH INTERFACE AREA REQUIRES AN EXTREMELY LIGHT TOUCH TO WORK EFFECTIVELY. THE LIGHTER THE FINGER PRESSURE, THE MORE LIKELY IT WILL SENSE THE TOUCH.

## 02 - SELECTING A WAVELENGTH

When the system is turned on, and the electronic passcode is properly input, you will be prompted with voice confirmation "Please select wavelength" and two flashing wavelength graphics to select the desired wavelength of choice. (Figure: 4.1) The Gemini EVO 810+980 soft tissue laser can operate in three wavelength modes: 810 nm alone, 980 nm alone, or Dual Wavelength. A wavelength mode must be selected before proceeding further, but may be changed at any time. (Figure: 4.2)

By selecting the desired wavelength, voice confirmation (if enabled) will sound as follows: "810 STBY", "980 STBY" "DUAL WAVELENGTH STBY"

When selecting a wavelength, the 3 wavelength options 810 / DUAL / 980 will be shown as they are the only necessary icons to be touched while selecting a wavelength.

## 03 - ACTIVATION PEDAL CONNECTION

- Connecting the activation pedal to your laser unit via Bluetooth for the first time is simple. (See Figures: 5.1–5.8)
- Install the provided 2 AA batteries into the Activation Pedal. (Figure: 5.1)
- Turn the laser unit ON. (Figure: 5.2)
- Enter the Passcode. (Figure: 5.3)
- Select the wavelength of choice. (Figure: 5.4)
- Press Active to initiate the Bluetooth connection between the laser unit and pedal. (Figure: 5.5)
- Depress and release the Activation Pedal once. Connection is done automatically. (Figure: 5.6)
- A Bluetooth indicator will appear on the display and Activation Pedal when it is properly connected and the Laser is in Active mode. (Figure: 5.7)
- The Activation Pedal comes with a protective shroud to prevent accidental laser activation. Please do not step on the protective shroud as it could result in accidental damage to the Activation Pedal. (Figure: 5.8)

## 04 - MANUAL POWER ADJUSTMENT

The Gemini EVO 810+980 soft tissue laser can output up to a maximum of 2.0 watts of average power. To adjust the power setting manually, touch the LEFT and RIGHT arrows on the Guided Touch Interface. Each touch of an arrow raises or lowers the power by 0.1 watts. Touching and holding an arrow will increase the speed in which the power setting is raised or lowered. Touch the ACTIVE button to put the laser in active mode. Depress the activation pedal to initiate the laser. (Figure: 6.1)

## CLINICAL TIP

Optimal results will be achieved by regulating the power and the speed that the operator moves the fiber optic tip. Tissue charring is an undesirable after-effect of too much power or of the fiber tip moving too slowly. Always use the least amount of power that is required to complete your procedure. The ideal tissue response will show little or no discoloration after treatment and will result in less collateral damage and faster healing.

Avoid penetrating or damaging the periosteum, and do not attempt to use the laser on alveolar bone. Because the laser energy is attracted to melanin and hemoglobin, power must be reduced when treating patients with darker pigmented soft tissue.

## 05 - LASER STBY AND ACTIVE MODES

The Active/STBY Guided Touch Interface selection serves a dual purpose. It activates (ACTIVE) and deactivates the laser (STBY). By default, the system powers up in STBY mode. The laser cannot be activated prior to selecting a wavelength. Each time the ACTIVE/STBY selection is touched, the system toggles between Active and STBY modes. There is an audio confirmation (unless voice confirmation is muted), and an icon for either "ACTIVE" or "STANDBY" near each wavelength indicator. The red aiming beam and tip illumination are visible only when the laser is in Active mode.

# CONTROLS, OPERATION & USE

When the system is in Active mode, touching any selection other than LEFT and RIGHT will return the system to the STBY mode. When the Activation Pedal is depressed in the Active mode, the outer indicator lines around each wavelength icon on the display flashes to provide a visual indication that the laser is firing. There is also an audio beep when laser is being fired. For safety purposes, a laser firing delay of 0.25 seconds was implemented in order to prevent accidental activation. (Figure: 6.2)

## 06 - PRESET PROCEDURE SETTINGS AND CUSTOMIZATION

Touch the PRESET selection to bring up all the preset procedures and categories on the display. A collection of pre-programmed procedures will be revealed on the display. Selecting the LEFT and RIGHT arrows will toggle between NON SURGICAL, SURGICAL, and PAIN RELIEF categories. Selecting the UP and DOWN arrows will toggle between procedures within each category. The corresponding power setting for each procedure is displayed on the Power Indicator when the procedure is highlighted. (Figure: 6.3, Option: 1)

A Simulated Screen showing the Gingivectomy Preset Procedure and dual wavelength mode selected. (Figure: 6.3)

The Gemini EVO 810+980 soft tissue laser is pre-programmed with 16 procedures listed under three categories: **NON SURGICAL**, **SURGICAL**, and **PAIN RELIEF**. Within each category are the most commonly used procedures with suggested power settings. Always use the minimum amount of power necessary to perform a particular procedure. Manual power adjustment may be necessary depending on patient and procedural needs.

NON SURGICAL	TISSUE CONTACT	810	DUAL	980
Decontamination	Yes	0.6w	0.5w	0.6w
Aphthous Ulcer	No	0.7w	0.6w	0.8w
Hemostatis	Yes	0.9w	0.8w	1.0w
Debridement	Yes	0.4w	0.3w	0.5w
SURGICAL	TISSUE CONTACT	810	DUAL	980
Incision / Excision	Yes	1.1w	1.0w	1.2w
Implant Recovery	Yes	1.3w	1.1w	1.5w
Tooth Exposure	Yes	0.8w	0.7w	0.9w
Operculectomy	Yes	1.2w	1.1w	1.4w
Gingivoplasty	Yes	0.8w	0.7w	0.9w
Gingivectomy	Yes	1.0w	0.9w	1.1w
Frenectomy	Yes	1.2w	1.1w	1.4w
Pulpotomy	Yes	0.8w	0.7w	1.0w
Troughing	Yes	0.9w	0.8w	1.0w
PAIN RELIEF	TISSUE CONTACT	810	DUAL	980
7 mm Tip	User Choice	0.3w	0.3w	0.3w
3 mm Tip	User Choice	0.3w	0.3w	0.3w
25 mm Tip	Spacer	1.0w	1.0w	1.0w

**Note:** Preset procedure settings are only a general recommendation from the manufacturer. They are not meant to replace the need for proper training or clinical judgement of the operator. The presettings are subject to changes through software updates and may therefore differ slightly from the settings indicated on this page. All power settings are shown in average power.

- Gemini EVO 810+980 soft tissue laser preset procedure settings can be customized. To save your own procedure setting, press PRESET (Figure: 6.3, Option: 2) once to bring up preset procedures on the display and navigate to the procedure you would like to customize.
- When a particular procedure is highlighted, press and hold PRESET ICON for 3 seconds.
- You will hear two audible beeps and the power indicator value and power bar will start flashing on the display. Use LEFT/RIGHT arrows to adjust the new average power to the desired setting.
- To save the setting, press and hold the PRESET ICON for 3 seconds again. You will hear two audible beeps when the setting has been saved.
- To reset all preset procedure settings to the factory default, press and hold the PRESET ICON for 5 seconds. You will hear three audible beeps when the settings are reset.
- Another way to customize the preset procedures is through the DASHBOARD.

**Hold 3 seconds = Customize / Save | Hold 5 seconds = Factory Reset**  
Power bar Indicator and power value will start flashing on the display. (Figure: 6.3, Option: 3)

## 07 - TIP ILLUMINATION

The Gemini EVO 810+980 soft tissue laser handpiece is equipped with a tip illumination light to provide better visibility of the surgical site during treatment. To toggle the light intensities between LOW, MEDIUM, HIGH, and OFF, touch the MENU icon and select TIP LIGHT on the Guided Touch Interface. Then use the arrows UP / DOWN to change LED intensities. LED will only stay on for 3 seconds as a preview when not in active mode. Tip Illumination icon shows Green color when this feature is enabled and Red when it is OFF.

Please note arrows UP / DOWN will appear and disappear according to the selected settings. As an example, if you select HIGH, the UP arrow will disappear indicating this is the highest setting available. The same behavior happens when you select OFF, in which the DOWN arrow will disappear. (Figure: 6.4)

## 08 - HAPTIC SENSE “HS” - ACTIVATION PEDAL & PBM

The Gemini EVO 810+980 soft tissue laser is equipped with Haptic Sense “HS” in order to add an additional dimension of feedback while depressing the activation pedal or during PBM procedures. This feature will provide the user a buzzing sensation to the foot or hand while the activation pedal is depressed. To toggle the Haptic Sense “HS” between activation pedal and handpiece (PBM procedures only), press the haptic icon. Blue icon indicates Haptic Sense is available for activation pedal and Amber icon for handpiece (PBM procedures only). To adjust intensities between LOW, MEDIUM, HIGH, and OFF, touch the arrows UP / DOWN to change the intensities. (Figure: 6.5)

# CONTROLS, OPERATION & USE

Blue - Activation Pedal (Figure: 6.5, Option: 1)

Amber - Handpiece for PBM procedures only (Figure: 6.5, Option: 2)

## 09 - SOUND

To change the sound level, touch the MENU and then the SOUND icon on the Guided Touch Interface. Adjust the sound level by touching the Up / Down arrows. To exit, touch the MENU icon to save your selection. The system remembers the last used sound setting when it is powered on. When UP arrow disappears, this indicates the volume is all the way to max and vice versa. Icon shows Red when OFF. (Figure: 6.6)

Voice Confirmation can be enabled and disabled by touching the voice confirmation selection on the Guided Touch Interface. Red icon shows disabled and Green enabled. (Figure: 6.6)

## 10 - AIMING LIGHT

To change the aiming light intensity, touch the MENU and then the AIMING icon on the Guided Touch Interface. Adjust the aiming light level by touching the Up / Down arrows. To exit, touch the MENU icon to save your selection. The system remembers the last used sound setting when it is powered on. When UP arrow disappears indicates the setting is all the way to max and vice versa. Icon shows Red when OFF. (Figure: 6.7)

## 11 - PHOTOBIMODULATION (PBM)

### WARNINGS & CAUTIONS

#### CAUTION:

- Do not connect or disconnect a PBM adapter while the Gemini EVO laser is turned on. Only connect or disconnect a PBM adapter when the Gemini EVO laser is inactive or standby mode.
- Do not use any harsh chemicals or abrasives to clean the glass optics within a PBM adapter. Doing so may damage the glass.
- Do not autoclave the 25 mm PBM adapter or spacers. Doing so will damage the components.
- The spacers are single-use only to avoid possible cross-contamination. They must be disposed of after use in a bio hazard medical waste Sharps container.
- Wavelength-appropriate eye protection must be worn at all times while using, and in proximity of, a PBM adapter while it is in use.

#### WARNING:

- The PBM adapters must only be used with a Gemini EVO laser. Do not attempt to use a PBM adapter with any other laser system or light source.
- Never look directly into a PBM adapter while the laser is active, even with safety eyewear on.
- Do not use the 25 mm PBM adapter without a spacer attached.

### PBM COMPONENTS

7 mm Adapter (Figure: 7.1, Option: 1) | 25 mm Spacer (Figure: 7.1, Option: 2) | 3 mm Adapter (Figure: 7.1, Option: 3) | 25 mm Adapter (Figure: 7.1, Option: 4)

### PBM ASSEMBLY

(Attachment threading procedure applies to all PBM tips equally)

1. Remove the dust covers
  - Dust cover / Front (Figure: 7.2, Option: 1)
  - Dust cover Back (Thread) (Figure: 7.2, Option: 2)
2. Screw the PBM adapter onto the end of the Gemini EVO laser handpiece until it is tight. (Figure: 7.2, Option: 3)
3. If using the 25 mm PBM adapter, screw a spacer onto the end of the 25 mm PBM adapter. (Figure: 7.2, Option: 4)

The PBM adapter is now ready to use. To remove the PBM adapter, unscrew it from the Gemini EVO laser handpiece and re-install dust covers when not in use.

### SELECTING AND ADJUSTING PBM PRESET

1. To enable Pain Relief, select PRESET on the Guided Touch Interface and navigate with the right arrow to the PAIN RELIEF category. (Figure: 7.3, Option: 1)
2. Use the UP / DOWN arrow select the desired PBM adapter. (Figure: 7.3, Option: 2)
3. Click ACTIVE to select treatment time in seconds. Display will flash 0.0 seconds. (Figure: 7.3, Option: 3)
4. Using the right arrow select the time in seconds followed by ACTIVE to enable timer. Press and hold right arrow advances timer faster in 10 seconds increments. (Figure: 7.3, Option: 4)

Laser unit is ready for PBM treatment. The timer counts down in seconds and stops automatically after treatment time is completed. If the activation pedal is released mid-treatment, the timer will pause and will resume when the pedal is depressed again.

### USE RECOMMENDATIONS

Affected muscles and/or joints have to be exposed to an adequate level of laser energy over a period of time to provide effective results. Some cases may require more than one laser treatment, or a series of treatments, before significant improvement is reported. Repeat the treatment as necessary and monitor the progress of the patient's condition throughout the treatment.

Diode laser wavelengths, especially 810 nm, are well absorbed in melanin in the skin, which can lead to greater heating of the target tissues in patients with darker skin types. Power and treatment time should be taken into account for patients of varying skin pigmentation. Refer to the Fitzpatrick Skin Type Scale for proper skin classification.

Pain relief preset procedure settings are programmed into the Gemini EVO laser for ease of use. Always use professional clinical judgment when selecting the laser settings for pain therapy.

Monitor the patient and adjust the power and/or treatment time as necessary to ensure both efficacy and patient comfort. The preset procedure setting is not meant to be a clinical recommendation in any way.

## CONTROLS, OPERATION & USE

When you are ready to begin treatment, hold the PBM adapter in contact with the target treatment area. The PBM adapter is designed to be held in a constant location for the duration of the treatment. If the desired treatment area is larger than the PBM adapter's spot size, move the adapter to a new location and start a new treatment only after the initial treatment time has elapsed.

### PBM ADVERSE EVENTS & CONTRAINDICATIONS

If patient discomfort or reddening of the skin in the treatment area occurs at any time during treatment, you can do the following:

- Defocus the laser energy by moving the adapter a few centimeters back from the skin
- Decrease the treatment time
- Stop treatment

If blistering of the skin occurs, or the patient feels a burning sensation, immediately stop treatment and rinse the area with cool water or place a cold pack on the affected area for at least 5 minutes. Afterwards, apply a burn ointment or spray. **DO NOT USE ICE.**

- Do not use over articles of clothing.
- Do not treat open wounds.
- Do not apply ointment, creams, lotions, or heating lotion patches at, or in close proximity to, the treatment area.
- Do not apply therapies prior to treatment that could change body temperature, such as ultrasound, ice/heat pack, electrical stimulation, or heating patches.
- Avoid treatment sites with tattoos.
- Different implant materials will respond differently to laser energy and heat; be aware of any implants and their location; avoid direct exposure to laser energy or heat at the site of the implant.
- Excessive fatty tissue is known to transmit heat without much attenuation, therefore increase distance or decrease treatment time.
- Muscle tissue closer to the skin surface may experience a higher absorption of heat; carefully monitor skin temperature and reduce treatment time as necessary.
- Patients with swelling and/or inflammation may be sensitive to heat; reduce treatment time as necessary to ensure comfort during treatment.
- Patients with tender or sensitive skin may be hypersensitive to heat; reduce treatment time as necessary to ensure comfort during treatment.
- Scar tissue has been associated with poor circulation and reduced cooling through heat transport by blood; reduce treatment time as necessary to avoid overheating.
- Do not treat directly over the site of any known primary malignant carcinoma or secondary metastasis except for palliative care with informed consent and oncologist permission.
- Do not treat pregnant women as the effects of photobiomodulation therapy on the fetus are unknown.

### PBM ADAPTER MAINTENANCE

The disposable spacers are supplied non-sterile by the manufacturer and should be wiped with isopropyl alcohol wipes by the operator prior to use. The spacers are intended for single-use only and should never be autoclaved or reused to prevent damage or cross-contamination.

The 25 mm PBM adapter is also provided non-sterile by the manufacturer and can be wiped as needed using isopropyl alcohol wipes. Do not submerge the 25 mm PBM adapter in any type of cleaning solution. **DO NOT AUTOCLAVE** the 25 mm PBM adapter.

Use the included cleaning cloth to gently wipe the glass optics of the 25 mm PBM as needed. Do not use any harsh chemicals or abrasives to clean the glass optics within the 25 mm PBM adapter. Doing so may damage the glass.

### PBM ADAPTER SPECIFICATION

	25 MM PBM ADAPTER	7 MM PBM ADAPTER	3 MM PBM ADAPTER
Dimensions	63 x 41 x 41 mm	100 x 24.7 x 14.9 mm	100 x 20.5 x 14.9 mm
Weight	76 grams	50 grams	50 grams
Spot Size	4.91 cm <sup>2</sup>	0.38 cm <sup>2</sup>	0.07 cm <sup>2</sup>
Power (under preset)	1.0 Watt Average	0.3 Watt Average	0.3 Watt Average
Power Density	204 mw/cm <sup>2</sup>	780 mw/cm <sup>2</sup>	4244 mw/cm <sup>2</sup>

### PBM DOSAGE TABLE

TIME IN SECONDS	25 mm PBM Dose (J/cm <sup>2</sup> )	7 mm PBM Dose (J/cm <sup>2</sup> )	3 mm PBM Dose (J/cm <sup>2</sup> )
03	0.6	2.3	12.7
06	1.2	4.7	25.5
09	1.8	7.0	38.2
12	2.4	9.4	50.9
15	3.1	11.7	63.7
18	3.7	14.0	76.4
21	4.3	16.4	89.1
24	4.9	18.7	101.9
27	5.5	21.1	114.6
30	6.1	23.4	127.3
33	6.7	25.7	140.1
36	7.3	28.1	152.8
39	8.0	30.4	165.5
42	8.6	32.8	178.2
45	9.2	35.1	191.0
48	9.8	37.4	203.7
51	10.4	39.8	216.4
54	11.0	42.1	229.2
57	11.6	44.5	241.9
60	12.2	46.8	254.6

# CONTROLS, OPERATION & USE

## 12 - DISPOSABLE TIP OPERATION

The disposable fiber tip is relatively flexible, but can be broken if bent at an angle that is too sharp. Use the provided Bending Tool to bend the tip to the desired angle. Do not bend the tip any more than the Bending Tool allows.

Protein debris from gingival tissue accumulates on the fiber tip during surgery, and the extreme heat that develops will deteriorate the optical efficiency. Fibers can fracture if a blackened area greater than 3-4 mm develops.

Replace the disposable, single-use fiber optic tip as necessary and for each new patient. The tips are provided in a sealed package. Each tip contains a pre-cleaved, pre-stripped piece of fiber. They are designed for single-use only and must be discarded after use. (Figures: 8.1–8.4)

DO NOT OVERBEND TIP / DO NOT BEND WITH FINGER

### DISPOSABLE TIPS

The Gemini's 5 mm single-use fiber tips are unique in that they come pre-initiated. That means there is black pigment added to the end of each fiber tip to help focus laser energy at the tip. All procedures that require the removal or cutting of soft tissue require an initiated tip. There are some procedures that call for an un-initiated tip, such as Aphthous Ulcer treatment where no tissue is being removed.

To ensure that the tip stays initiated when wiping the tip with isopropyl alcohol before a procedure, activate and fire the laser at 1 watt of average power for 1-2 seconds prior to wiping the tip. This action will ensure that the pre-initiation does not wipe off during the cleaning process.

When a procedure calls for an uninitiated tip, simply rub off the pigment at the end of the Fiber Tip with gauze and isopropyl alcohol. This pigment removal must occur before firing the laser.

The Gemini's 7mm single-use fiber tips come uninitiated. Laser procedures that do not remove tissue, such as decontamination or the treatment of aphthous ulcers, do not require the laser tip to be initiated. One way to initiate the tip is to rub the tip on articulating film while firing the laser on a low power setting. (Figure: 8.5)

**IMPORTANT NOTE:** Not all soft-tissue procedures require an initiated—or darkened—tip and contact with tissue. Those procedures that do not require tissue contact will use a fiber tip that is NOT initiated, because to be effective in non-contact mode, laser energy must flow unimpeded from the tip into the target tissues. Follow the procedure above to uninitiate a fiber tip.

### 5MM TIPS (PRE-INITIATED)

Surgical procedures such as Incision/Excision, Implant Recovery, Tooth Exposure, Operculectomy, Gingivoplasty, Gingivectomy, Frenectomy, and Troughing are some of the procedures recommended with a 5 mm tip. (Figure: 8.6)

### 7 MM TIPS (UN-INITIATED)

Decontamination and Aphthous Ulcer are some of the procedures recommended with a 7 mm tip. (Figure: 8.7)

## 13 - BATTERY AND BATTERY LEVEL INDICATIONS

The Gemini EVO 810+980 soft tissue laser is equipped with a powerful Lithium-Ion battery capable of delivering a full day of laser usage. Simply connect the provided power supply to the rear of the unit and charging will start immediately.

- It is recommended to fully charge the laser unit prior to initial use after unpacking.
- The battery level indicator is located at the upper right corner of the display and shows battery percentage remaining.

 100%	<b>Usage Time:</b>	2.0 Hours of consecutive max power of 2.0W
 75%	<b>Usage Time:</b>	1.5 Hours of consecutive max power of 2.0W
 50%	<b>Usage Time:</b>	1.0 Hour of consecutive max power of 2.0W
 25%	<b>Usage Time:</b>	30 min of consecutive max power of 2.0W
 0%	Minimum of 60 minutes of charging is required before first use	

 "Please connect power supply" reminder.

- Standby time refers to when the unit is off and not being used.
- Usage time refers to constant use without interruption.
- To preserve battery life, the laser unit goes into "Inactivity" mode within 10 minutes of inactivity.
- The Lithium-Ion battery has a typical lifespan of 2 years, at which time it is advised that the battery be replaced.

## 14 - POWER SUPPLY

Use only the provided 13V, 4A AC/DC power supply for charging the system battery and as an alternate laser power source. During initial setup use the AC/DC power supply for one hour to fully charge the battery.

Plug the power supply into an AC outlet and connect to the corresponding connector on the rear of the laser Unit. Only use the power supply provided with the system.

**WARNING:** To avoid the risk of electric shock, this equipment must only be connected to a supply main with a protective earthing conductor.

**ONLY USE THE 18V POWER SUPPLY WITH THE Gemini EVO 810+980 DIODE LASER. OTHER POWER SUPPLIES INCLUDING THE POWER SUPPLY FROM OTHER Gemini Laser PRODUCTS CAN CAUSE DAMAGE TO YOUR Gemini EVO Laser UNIT.**

## 15 - FIBER WRAPPING

A fiber wrapping system was built within the laser unit in order to provide a safe and convenient way to manage and store the optical fiber system. To store the fiber properly, always wrap in a clockwise direction to protect and store the fiber optic cable when not in use. (See Figures: 9.1–9.4)

The fiber optic cable conducts laser energy from the laser diodes to the target tissues. These fibers are made of a thin glass silica. Note that there are potential hazards when inserting, steeply bending, or improperly securing the fiber optic tips to the handpiece. Failure to follow these recommendations may lead to damage of the fiber or delivery system, and/or harm to the patient, staff, or laser operator.

# CONTROLS, OPERATION & USE

**CAUTION: AVOID DAMAGING THE FIBER.** Do not wrap the fiber in a counter-clockwise direction. Doing so will possibly damage the optical fiber, preventing the use of the laser.

## 16 - HANDPIECE MAGNET

The Gemini EVO 810+980 soft tissue laser is designed with a strong magnet that will secure the surgical handpiece in place when the laser is not in use. Gently place the handpiece behind the transparent display or over the "rounded pad" of the laser unit and the magnet will hold the handpiece in place. (Figure: 10.1)

- Rounded rubberized pad (Figure: 10.1, Option: 1)
- Magnets (Figure: 10.1, Option: 2)

## 17 - OPERATING MODE

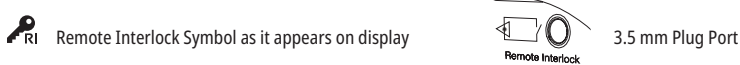
The Gemini EVO 810+980 soft tissue laser will only deliver energy in Pulsed "temporal emission mode" and is optimized to provide the operator ideal control of target tissue temperatures and efficiency of energy delivered. The pulse width is fixed and not user adjustable. The operator will only need to adjust the laser wavelength and average power.

## 18 - EMERGENCY STOP

The Gemini EVO 810+980 soft tissue laser can be immediately deactivated in any mode, at any time, and in any power setting by pressing the red STOP button located in the front left of the system. (Figure: 1.3, Option: 6)

## 19 - REMOTE INTERLOCK (SWITCH NOT INCLUDED)

The Gemini EVO 810+980 soft tissue laser provides a remote interlock feature that enables a clinician to establish a dedicated laser treatment room with a remote interlock connector. A switch on the entrance door is attached and is electronically wired into the Laser Unit via 3.5 mm plug jack. When the door to the room is opened, the connector/switch provides an electrically open circuit that deactivates laser emissions. To use the remote interlock feature, an interlock connector/switch and cable must be purchased. Contact the manufacturer for assistance. (Figure: 1.4, Option: 6)

## 20 - ADD NEW ACTIVATION PEDAL - BLUETOOTH CONNECTION

In certain instances, your office may need to have an additional Activation Pedal with the same Gemini EVO device. You can order an additional Activation Pedal and follow the instructions below to properly pair the new pedal with your Gemini EVO unit. For a successful pairing, the sequence below must be followed:

1. Depress the Pedal (Figure: 11.1)
  - If batteries are installed, remove both of them. Keep the pedal depressed until completion of STEP 3.
2. Install Batteries (Figure: 11.2)
  - Keep the pedal depressed. Install the 2 batteries and keep the pedal depressed for 15 seconds.
3. Flashing Icons (Figure: 11.3)
  - Activation Pedal icons will flash simultaneously along with the haptic feedback indicating pairing mode sequence has started. You can now release the pedal. (Figure 1.5, Options: 7-10)
4. Pairing
  - Press and hold the HOME icon for 5 seconds.
  - "Bluetooth Pairing Enabled" (Figure: 11.4)
  - Activation Pedal automatically connects to your Gemini EVO device.
  - "Bluetooth Pairing Successful" (Figure: 11.5)

## 21 - TRANSPARENT ELECTROLUMINESCENT DISPLAY

The Gemini EVO 810+980 soft tissue laser is designed with a unique transparent electroluminescent display that can provide extremely high resolution viewing angles from up to 160 degrees of field of view. The transparent display was specially designed with over 80% transparency and an arch of 15 degrees for optimal viewing angle from any direction. The light is generated by a thin film, less than 2 microns thick, of specially designed electroluminescent phosphor.

- **Do Not Grab unit by display** (Figure: 12.1)
- **Grab unit by the base** (Figure: 12.2)

## DISPLAY SPECIFICATIONS

Display type:	Electroluminescent
Transparency:	80%
Brightness:	300 cd/m <sup>2</sup>
Color:	Broadband Yellow
Peak Wavelength:	582 nm
Voltage:	~195V AC
Response time:	1.8 ms
Glass type:	Soda lime
Glass thickness:	4.1 mm fused
Thin Film Thickness:	~2 micron

# MOBILE APP & DASHBOARD

## 22 - ENABLING WI-FI CONNECTIVITY VIA APP

The Gemini EVO 810+980 soft tissue laser enables you to bridge your location's existing Wi-Fi network directly to the Gemini EVO unit. This provisioning enables your Gemini EVO unit to receive internet connection, allowing the user to receive important performance updates, technical support, track procedures, and a multitude of other features.

To properly enable Wi-Fi connectivity, please follow the steps below:

# MOBILE APP & DASHBOARD

## 1. DOWNLOAD THE MOBILE APP

The Gemini EVO Laser App is available for iOS and Android devices. On your mobile app store, search for Gemini EVO laser. (Figure: 2.1)

## 2. REGISTERING AND AUTHENTICATING YOUR DEVICE

Once you have the App installed, select "I need an account" and follow the easy step-by-step instructions to register your laser. Registering your unit is an important step to allow your Gemini EVO laser to receive internet connection. You will receive an authentication code via email. (Figure: 13.1)

## 3. SCAN YOUR LASER

Your Gemini EVO unit contains a unique QR code located on the bottom of your unit or activation pedal. Point your phone camera to the QR code and the app will scan the laser unit. You can nickname your Gemini EVO laser any name you like and click SAVE. (Figure: 13.2)

## 4. ENABLING WI-FI

After saving your device name above, follow the simple steps on the App in order to enable the laser unit to start talking to your local Wi-Fi network. (Figure: 13.3–13.4)

1. Turn the Gemini EVO unit ON
2. Enter passcode
3. Select wavelength (any wavelength)
4. Press and hold the Wi-Fi icon for 3 seconds. "Wi-Fi Setup in Progress" voice confirmation is heard. Move to next step to select a local Wi-Fi.

## 5. SELECTING A WI-FI NETWORK

A list of Wi-Fi networks will be displayed. Please select the Wi-Fi network associated to your office location and enter appropriate password. Please note the Gemini EVO Laser Wi-Fi is compatible with 2.4 GHz networks only. If you have a secure firewall or Anti-Virus software, you may need to contact your network administrator in case there are difficulties connecting to your local Wi-Fi network. (Figure: 13.5)

## 6. ESTABLISHING A WI-FI CONNECTION

After selecting the appropriate Wi-Fi network and entering the password, the Gemini EVO unit will establish a secure connection with your local Wi-Fi. The Wi-Fi connection between your local Wi-Fi network to the Gemini EVO unit can take up to 2 minutes to complete. The Gemini EVO unit displays a progress bar on the electroluminescent display. Upon connection the App displays the main page shown below. (Figure: 13.6)

### SEE FIGURE 13.7

- |  |   |
|--|---|
| 01 - Gemini EVO laser name given during registration                                     | 07 - Edit name registered   |
| 02 - Wi-Fi connectivity status: Green = Connected to Wi-Fi, Red = Not connected to Wi-Fi | 08 - Overall Laser unit information: Green = Receiving status, Red = Not receiving status |
| 03 - Latest available software update  | 09 - Extended warranty availability and status  |
| 04 - Add / Remove Gemini EVO devices   | 10 - Usage statistics   |
| 05 - Videos, Instructions, Techniques  | 11 - Up to date User's Manual   |
| 06 - Log in and out different users  | 12 - Account info / Change personal information   |

Similar to identical user interface for Android users. Layout is subjected to change based on future updates on all platforms.

### UPDATES (Figure: 14.1)

With the iOS and Android App, you can perform automatic updates directly to your Gemini EVO laser. Automatic updates are extremely important as it enables your Gemini EVO laser to utilize the latest and greatest improvements.

### DEVICES (Figure: 14.2)

This page will enable you to add or remove a Gemini EVO device from your registered account. You can have multiple Gemini EVO devices registered with one account. Units shown in green are currently online. Units shown in red are currently offline.

### VIDEOS (Figure: 14.3)

The videos tab will show you several of the procedures that can be performed with the Gemini EVO laser. Additionally, we will upload the latest techniques and customer tips as reference.

### STATUS (Figure: 14.4)

The status tab will show several important status conditions such as the health of your battery, the strength of your Wi-Fi connection, and the ability to add and remove another Gemini EVO laser from your account. The Status page is the overall health of your Gemini EVO device.

### WARRANTY (Figure: 14.5)

With the iOS and Android apps, you can purchase extended warranty if qualified. The Gemini EVO laser comes with 2 years limited factory warranty. You can extend your factory warranty to an additional 24 or 36 months. Warranty starts from shipment day. Prices shown are subjective to change.

### STATISTICS (Figure: 14.6)

With the Gemini EVO app, you can track how many procedures have been performed by category, see which wavelength mode is used the most, as well as the overall laser usage time of this Gemini EVO device.

### USER MANUAL (Figure: 14.7)

With the iOS and Android apps, you can have access to the user's manual at any time. The user's manual will always contain the latest update, enabling you to always have access to the latest documentation.

### ACCOUNT (Figure: 14.8)

The account page allows you to change your registered name, phone number, and product nickname. This is an important feature in case the Gemini EVO device is exchanged with another office.

# MOBILE APP & DASHBOARD

## WEB INTEGRATION VIA THE DASHBOARD

Once connected to Wi-Fi, the Gemini EVO 810+980 soft tissue laser will share data with the DASHBOARD, which will enable you to visualize several parameters of your laser. Login to [dashboard.geminievo.com](https://dashboard.geminievo.com) and use the same login credentials created within the App for iOS and Android devices. (Figure: 15.1)

## GENERAL USE

The Dashboard main page enables you overall visibility of your Gemini EVO laser unit. Some of the information displayed on the dashboard includes:

- Graphical representation of manual or preset procedures performed
- The most used power settings
- Total number of procedures performed by day / week or a specific time range selected by user
- User's total laser usage
- Wavelength usage chart showing the most used wavelengths
- Top preset procedures performed
- Top preset procedures performed by category (Non Surgical / Surgical / Pain Relief)
- Return on investment visibility
- Also access additional support resources

The Dashboard is being constantly improved. Some of the features listed above might be different and updated/improved since product inception. Our goal is to always improve the system based on customer feedback. If you have an improvement suggestion, please email to [feedback@azenamedical.com](mailto:feedback@azenamedical.com) and we will do our best to analyze and implement it into our next update.

## DASHBOARD MENU

The main menu of the dashboard contains several links that facilitates the navigation of all the features listed, including the PBM Calculator which lets you visualize PBM treatments based on pain level, tissue color, or time. (Figure: 15.2)

## DASHBOARD

The main dashboard page contains all the up to date global information of your Gemini EVO device.

## CHAT SUPPORT

We have incorporated a very simple and user friendly way to get technical support. Chat and video are available during our normal business hours for your convenience.

## YOUR DEVICES

With one or several Gemini EVO devices, you can monitor and adjust several of your Gemini EVO laser features such as battery health under vitals or adjust preset procedures under the customize page. Some additional features also available on the dashboard are listed below:

### **Preset vs Manual**

This feature shows graphically the procedures performed manually vs with presets. This is a good way to visualize which procedure method the user is more comfortable using.

### **Top Procedures**

This feature shows the top procedures performed by name and by category. This is a good way to visualize which procedure is most performed by the user.

### **Total Procedures**

This feature shows the total amount of procedures performed in a specific period of time.

### **Return on Investment**

Upon adding the cost of each procedure based on user's region, the ROI feature will automatically calculate and display the return made on each procedure.

### **Wavelength Usage - Global**

This feature will display graphically which wavelength the user has used the most.

## DASHBOARD - MAIN PAGE - GLOBAL DATA

Once connected to your Gemini EVO laser, the Dashboard main page will display status information on your laser device. The Dashboard page will always display your overall data, even if you have multiple Gemini EVO devices associated with the same account. For instance, under Total Procedures, it will list all procedures performed with all your Gemini EVO devices. If you have one Gemini EVO device, it will list information about only that device. (Figure: 15.3)

## DASHBOARD - VITALS - UNIT SPECIFIC

The vitals page will display specific information about your Gemini EVO device by unique serial number. This type of information is helpful while troubleshooting or visualizing a feature when the user has a functionality issue. Some of the information that will be displayed is Top Procedures by category, most used wavelength, components status, and battery health. (Figure: 15.4)

## DASHBOARD - CUSTOMIZE

The customization page enables you to customize settings on your Gemini EVO device. In addition to changing several of the settings, you can purchase an extended warranty when available and perform software updates when available. (Figure: 15.5)

### **Customizing Preset Procedures**

Customizing a factory preset power setting to any other power settings is simple. Change the power by pressing the + or - icon and click save. Your unit will beep when the new preset is saved. To reset power to factory settings, simply click "Factory Reset" and confirm your selection.

### **Changing your Gemini EVO device name**

Your Gemini EVO laser name can be changed to any name you like. Select the edit icon under Laser Units registered and change the name. Select save to register the change.

### **Software Updates**

We will periodically release software updates. The dashboard will display if there's a software update available. Simply click update and your unit will automatically update to the latest release.

### **Customize ROI**

The ROI tool is a simple way for you to check how much revenue the Gemini EVO laser has generated. Add the cost of each procedure based on your region, and the dashboard will automatically calculate based on the amount of procedures you have performed. The ROI will be displayed in the main dashboard page.

### **Warranty**

Extended warranty (24 or 36 months) can be purchased via the dashboard.

## **CLEANING AND STERILIZATION PROCEDURES**

### **GUIDELINES**

The Gemini EVO 810+980 soft tissue laser is not supplied in sterile condition, nor must it be sterilized before use with the exception of the handpiece. The following cleaning and sterilization procedures are recommended before the initial use, and after each subsequent use:

1. The disposable fiber tips are supplied non-sterile by the manufacturer and should be wiped with isopropyl alcohol wipes by the operator prior to use. The tips are to be discarded in an infectious waste container (Sharps) after each use. There is no re-use or re-processing procedure indicated for the disposable fiber tips.
2. The aluminum handpiece, 3 mm, and 7 mm PBM adapters are also provided non-sterile by the manufacturer and should be cleaned and sterilized prior to initial use and after each use following these instructions:

### **HANDPIECE CLEANING**

**Warning: The Gemini EVO 810+980 soft tissue laser or any of its components cannot be cleaned with an automated cleaning process.**

The cleaning process is intended to remove blood, protein, and other potential contaminants from the surfaces and crevices of reusable accessories. This process may also reduce the quantity of particles, microorganisms, and pathogens present. Cleaning must be performed within a maximum of 1 hour after the procedure and always prior to sterilization:

1. After use, carefully remove the disposable fiber tip from the handpiece and dispose of in an infectious waste container (Sharps).
2. Clean the handpiece and attached fiber cable by using one CaviWipes® towelette, or equivalent product, to completely pre-clean exposed areas of all gross debris. Be sure to wipe the threaded area where the disposable tip attaches. Same procedure applies for 3 mm and 7 mm PBM adapters. PBM adapters must be removed from handpiece prior to cleaning.
3. Remove the handpiece shell following the instructions on the next page. Use a fresh towelette to pre-clean the newly exposed handpiece shell thread and the end cap (starting at the O-ring and wiping towards the fiber cable).
4. Use a new towelette to thoroughly wet all pre-cleaned areas, keeping all areas wet for 2 minutes at room temperature (68° F/20° C). Repeated use of towelettes may be required to ensure that the surfaces remain visibly wet.
5. Visually inspect the handpiece to ensure there is no visible debris remaining. If necessary, continue wiping with CaviWipes® until all visible debris is removed.
6. Wipe all exposed areas of the handpiece with isopropyl alcohol wipes to remove any residue left by the CaviWipes®.

### **HANDPIECE STEAM STERILIZATION**

The steam sterilization process is intended to destroy infectious microorganisms and pathogens. Always perform the sterilization procedure immediately after cleaning and prior to use, and only use FDA-cleared (USA) or CE-marked (Europe) sterilization accessories such as sterilization pouches and autoclave trays. (Figure: 17.1)

1. Place the handpiece shell, 3 mm and/or 7 mm PBM adapters in a separate single-wrap, self-seal autoclave pouch.
2. Place on an autoclave tray with paper side down; do not stack other instruments on top of the pouch.
3. Place the tray inside the autoclave chamber and set the cycle to 135° C (275° F) for a minimum of 10 minutes, with a dry time of 30 minutes.
4. Once the cycle is completed, remove the tray and let the sterilized item cool and dry. The handpiece, 3 mm and/or 7 mm PBM adapters must remain in the sterilization pouch until used in order to maintain sterility.
5. Visually inspect handpiece shell or 3mm/7mm PBM adapter to ensure product is not degraded. Below are the criteria for the degradation for the respective items:

A visual and mechanical inspection of the PBM adapters and aluminum handpiece after each sterilization should be conducted to ensure adapters have not degraded and lost performance.

Unacceptable deterioration includes cracked glass (Figure: 16.1), delamination of anodized material (Figure: 16.2), a non-circular spot when checking the aiming light on a flat surface (Figure: 16.3) and not being able to fully thread onto the handpiece. In the event the adapters have cracked glass or a non-circular aiming light spot, please send the adapters back to the manufacturer for review.

6. Remove / reassemble the handpiece shell or 3mm/7mm PBM adapter following the instructions below.

- 7mm PBM Tip turn counter clockwise to remove it. (Figure: 17.1, Option: 1)
- 3mm PBM Tip turn counter clockwise to remove it. (Figure: 17.1, Option: 2)
- Turn aluminum handpiece shell counter-clockwise to be removed. (Figure: 17.1, Option: 3)
- Remove aluminum handpiece shell for cleaning and sterilization. (Figure: 17.1, Option: 4)
- Extreme care is needed not to accidentally damage LED lens. (Figure: 17.1, Option: 5)
- Clean from these threads down the fiber optic cable. (Figure: 17.1, Option: 6)
- To reassemble, slide the handpiece shell onto the handpiece body and turn clockwise to tighten.

**NOTE:** The exterior of the laser unit is not routinely contaminated by procedures. The Guided Touch Interface and Electroluminescent display should be covered with a protective clear adhesive barrier film, replaceable after each patient. If the exterior of the laser unit should become contaminated, it should be wiped down with CaviWipes®, or equivalent product, then re-covered with a new protective plastic cover. We recommend wringing out any cleaning wipes before use to avoid dripping liquid on the laser unit.

**DO NOT** spray any disinfectant directly on the laser unit, because it could damage the transparent electroluminescent display.

**DO NOT** use abrasive materials to clean the laser or the display.

Never point the Laser tip directly at the face, eyes, or skin of anyone while emitting energy.

## **PROCEDURAL RECOMMENDATIONS**

### **GUIDELINES**

The following procedure guidelines are provided as a guide only and have been developed based on information provided by experienced laser users and educators. Always review the patient's history to evaluate possible contra-indication for use of local anesthesia or other complications.

# PROCEDURAL RECOMMENDATIONS

All clinical procedures performed with the Gemini EVO 810+980 soft tissue laser must be subjected to the same clinical judgment and care as with traditional techniques and instruments. Patient risk must always be considered and fully understood before clinical treatment. The clinician must completely understand the patient's medical history prior to treatment.

## INDICATIONS FOR USE

The Gemini EVO 810+980 soft tissue laser is intended for the incision, excision, ablation, vaporization, hemostasis, and treatment of oral soft tissue.

- Excisional and incisional biopsies
- Exposure of unerupted teeth
- Fibroma removal
- Gingivoplasty
- Hemostasis and coagulation
- Incision and drainage of abscess
- Operculectomy
- Pulpotomy
- Reduction of gingival hypertrophy
- Soft tissue crown lengthening
- Vestibuloplasty
- Laser Soft Tissue Curettage
- Tissue retraction
- Frenectomy and Frenotomy
- Gingival troughing for crown impressions
- Gingivectomy
- Gingival incision and excision
- Implant recovery
- Leukoplakia
- Oral papillectomies
- Pulpotomy as an adjunct to root canal therapy
- Reduction of bacterial level (decontamination) and inflammation
- Treatment of aphthous ulcers
- Lesion (tumor) removal
- Treatment of canker sores and herpetic and aphthous ulcers of the oral mucosa
- Removal of Diseased, Infected, Inflamed and necrotic soft tissue within the periodontal pocket
- Removal of highly inflamed edematous tissue affected by bacteria penetration of the pocket lining and junctional epithelium
- Sulcular debridement (removal of necrotic, diseased, or

inflamed soft tissue in the periodontal pocket to improve clinical indices including gingival index, gingival bleeding index, probe depth, attachment loss, and tooth mobility)

- Pain therapy. Topical heating for the purpose of elevating tissue temperature for a temporary relief of minor muscle and joint pain and stiffness, minor arthritis pain, or muscle spasm, minor sprains and strains, and minor muscular back pain, the temporary increase in local blood circulation; the temporary relaxation of muscle.

**All procedures listed in this manual are safe if performed by a licensed, trained professional. The potential side effects to the patient can include swelling, inflammation, redness of the skin, scarring, tissue pigment changes, and infection after treatment. All of these conditions can be reduced by cautiously following the appropriate aftercare or post-operative care instructions.**

# FACILITY & ENVIRONMENTAL CONSIDERATIONS

## GUIDELINES

In addition to receiving proper training in the use of soft-tissue dental lasers, users should be familiar and experienced with these procedures using electrosurgical devices or traditional instruments before performing them on patients with the Gemini EVO 810+980 soft tissue laser. Non-experienced users should seek appropriate training guidance before attempting clinical treatments with the Laser system. In order to insure the safe use of the Gemini EVO 810+980 soft tissue laser in your facility, please check to make sure that the proposed location is compatible with the specifications listed below.

## POWER REQUIREMENTS

External AC/DC Power Supply - Use only the provided Gemini EVO laser power supply. Every Gemini EVO laser power supply shows the corresponding label below. DO NOT use any other power supply. (Figure: 18.1)

Input Power: 100-240V; 50-60Hz, 1.5A  
Output Power: 18V, 65W

## HEATING AND VENTILATION

Operating environmental conditions to be within 10° - 40°C (50° - 104°F), and 95% relative humidity or less. Transportation and storage environmental conditions to be within 0° - 40°C (32° - 104°F), and relative humidity of 95% or less. Atmospheric pressure to be within 70kPa - 106kPa in operating, transportation and storage conditions.

## COMBUSTIBLE CHEMICALS AND GASES

All gases that are combustible or support combustion and are used in the operatory area where the Gemini EVO 810+980 soft tissue laser is being operated must be turned off during the procedure. Cleaning supplies or other flammable chemical compounds should be stored in an area away from the surgical site in order to avoid possible combustion. Do not use in the presence of supplemental therapeutic oxygen supplies for patients with respiratory or related diseases.

## PLUME EVACUATION

Plume evacuation should be addressed when vaporizing tissues. A high volume vacuum system should be used and 0.1 micron or less high filtration masks that are suitable for virus and bacterial control should be worn by clinicians.

## OPERATORY ACCESS DURING LASER USE

Access to the treatment area should be restricted while the lasers are in use. A sign indicating "LASER IN USE" should be placed in a designated area adjacent to the treatment area entry location.

**CAUTION: Laser fume and/or plume may contain viable tissue particulates**

# GENERAL SAFETY CONSIDERATIONS

## GUIDELINES

Safe use of the Gemini EVO 810+980 soft tissue laser is the responsibility of the entire dental team including the doctor, any system operators, and the dental office safety officer. In order to properly assess the favorable conditions of treatment, below is a pre-treatment checklist to help ensure treatment to your patient is safe:

- Ask the patient about allergy to local or topical anesthetics.
- Make sure the Laser Warning sign posted in the operating area.
- Make sure the patient and operator(s) are all wearing laser protective eyewear specific to Gemini EVO laser.
- Have the patient fill out an informed consent form for laser treatment. Form templates are typically available from your laser training provider.

# GENERAL SAFETY CONSIDERATIONS

- If performing a non-surgical procedure, use an un-initiated fiber tip.
- If performing a surgical procedure, use an initiated fiber tip.

Adjust the laser power settings as needed to fit the clinical circumstances of the case. The preset procedure settings built into Gemini EVO laser are simply a manufacturer’s recommendation. Optimal power level may vary case by case

## CONTRA-INDICATIONS

Exercise caution for general medical conditions that might contra-indicate a local procedure. Such conditions may include allergy to local or topical anesthetics, heart disease, lung disease, bleeding disorders, and immune system deficiency, or any medical conditions or medications that may contra-indicate use of certain light/laser type sources associated with this device. Medical clearance from patient’s physician is advisable when doubt exists regarding treatment.

The Gemini EVO 810+980 soft tissue laser is not indicated for hard tissue procedures. The Laser is attracted to melanin, hemoglobin, and to some extent, water. Avoid prolonged exposure of the energy when working in and around the cervical areas of the tooth. Due to the thin layer of enamel in this area, energy may be absorbed by the hemoglobin in the pulp and pulpal hyperemia may occur. Extended exposure to such energy could cause patient discomfort and even lead to possible pulpal necrosis.

## EYE AND SKIN PROTECTION

While the Gemini EVO 810+980 soft tissue laser is in use, doctors, system operators, auxiliary staff, patients, and anyone in the operatory must wear the appropriate safety eyewear that has been designed for use with the 800-plus nm wavelengths associated with lasers. Eye protection must conform to Specification DIN EN207 Annex II of the Directive 89/686/EEC with optical density of OD+5 for the wavelength range of 800nm-1000nm.

Nominal Ocular Hazard Distance (NOHD) is the distance from the source of laser emission to the point where it no longer exceeds its Maximum Permissible Exposure (MPE – highest level of laser radiation to which a person may be exposed without hazardous effects or adverse biological changes in the eyes or skin). The Nominal Hazard Zone (NHZ) is the space within which the level of direct, reflected, or scattered radiation during normal operation exceeds the appropriate MPEs. The outer limit of the NHZ is equal to the NOHD. The NOHD for persons NOT wearing recommended safety glasses is shown in Table 1 below.

TABLE 1: NOHD (INCHES/CM)

RADIATION SOURCE	MPE $\mu\text{j}/\text{cm}^2$	Divergence Angle	Without Eye Protection
FIBER OPTIC TIP	1.50	14.5° (+/- 1°)	77 in / 196cm
PBM 25 mm	1.50	4.5° (+/- 1°)	244 in / 620cm
PBM 7 mm	1.50	68° (+/- 1°)	16 in / 42cm
PBM 3 mm	1.50	68° (+/- 1°)	16 in / 42cm

**Never point the laser tip directly at the face, eyes, or skin of anyone while emitting energy.**

## EMERGENCY SHUTDOWN OPTIONS:

Perform any of these actions to terminate laser emissions in the event of a real or perceived emergency:

Press the emergency “STOP” Button. (Figure: 19.1)

Press the “ON/OFF” button. (Figure: 19.2)

Remote interlock open circuit deactivates the Laser. (Remote Interlock switch provided by request) (Figure: 19.3)

Touch the ACTIVE/STBY Guided Touch Interface selection (Figure: 19.4)

Release your foot from the Activation Pedal. (Figure: 19.5)

# SYSTEM SPECIFICATIONS

## GEMINI EVO 810+980 SOFT TISSUE LASER

<b>Dimensions of Foot Pedal:</b>	6.4" (L) x 4.2" (W) x 5.0" (H) - 16.5cm(L) x 10.8cm (W) 13.0m (H)
<b>Weight Laser Unit:</b>	2.8 lbs - 1.2Kg      Weight Activation Pedal: 0.6 lbs - 0.2Kg
<b>Laser classification:</b>	Class IV laser device
<b>Delivery system:</b>	Optical Fiber
<b>Wavelength:</b>	810nm or 980nm $\pm$ 10nm Dual Wavelength $\pm$ 10nm (50% @ 810nm / 50% @ 980nm)
<b>Maximum average power:</b>	810nm @ 2.0 Watts $\pm$ 20% 980nm @ 2.0 Watts $\pm$ 20% Dual Wavelength @ 2.0Watts $\pm$ 20%      Peak: Up to 100W
<b>Aiming beam wavelength:</b>	650 $\pm$ 10nm
<b>Aiming beam power:</b>	2mW max
<b>Beam divergence:</b>	254 mrad
<b>Power range:</b>	0.1 Watt to 2.0 Watts Average
<b>Pulse frequency:</b>	20 - 800Hz
<b>Pulse width:</b>	0.05ms
<b>Duty cycle:</b>	0.1% - 4%

# SYSTEM SPECIFICATIONS

<b>Voice confirmation:</b>	YES
<b>Power requirements:</b>	Input: 100-240 VAC @ 50 to 60 Hz - 1.5A Output: 1.8V DC, 3.6A, 65W
<b>Battery:</b>	14.4V Rechargeable Lithium Ion
<b>Wireless frequency:</b>	Bluetooth at 2.4 GHz Wi-Fi: 2.4 GHz / WPA2, WPA-Enterprise
<b>Maximum Operating Altitude:</b>	5,000 meters or 16,404 feet
<b>The Gemini EVO 810+980 SOFT TISSUE LASER complies with the following: IEC 60825-1 / EN/S 60601-1 IEC 60601-1-2 IEC 60601-2-22 21 CFR 1040.10 and 1040.11 FCC parts 15 and 18 (47 CFR)</b>	

## SERVICE AND TROUBLESHOOTING

### CALIBRATION

Re-calibration is recommended every 12 months in order to assure accuracy of optical output power. The Gemini EVO 810+980 soft tissue laser may be returned to the manufacturer for re-calibration, which you can arrange by contacting your distributor. Certain government or corporate entities may require calibration certificates which can also be provided by the manufacturer.

### ADVERSE EFFECTS

If used properly, there are no known adverse effects of using the Gemini EVO 810+980 soft tissue laser. Please thoroughly read and understand all Warnings, Precautions, and Contra-indications in this manual prior to use. In the event the laser malfunctions due to exposure to certain environment conditions, magnetic fields, external electrical influences, electrostatic discharge, pressure or variations in pressure, accelerations, and any potential thermal ignition sources, discontinue use and follow the instructions in the service and troubleshooting section of this manual. Additional measures may be necessary such as reorienting or relocating the device.

No separate equipment is recommended to be used to assess the favorable conditions which are acceptable for treatment or to assess the unfavorable conditions which would render a treatment unacceptable or hazardous.

Maximum LASER OUTPUT of laser radiation with the magnitudes of cumulative measurement uncertainty and any expected increase in the measured quantities after manufacture is stated as the standard uncertainty of measurement.

### WIRELESS INTERFERENCE

This equipment has been tested and found to comply with the limits for Class B Digital Device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. To determine if this equipment causes harmful interference to radio or television reception, turn the equipment off and on.

### ALL OTHER CONDITIONS

In the event that the Gemini EVO 810+980 soft tissue laser fails to operate correctly, and your distributor representative is unable to help, the system will need to be returned to the manufacturer for repair. There are no user-repairable parts available for the device. It is recommended that the system be returned in its original shipping box. If not available, one can be requested at the time you discuss your service incident with your distributor representative.

## TROUBLESHOOTING GUIDE

Why is it difficult to see the red aiming light or why is it barely visible?		My activation pedal is not working	
<b>CAUSE:</b>	<ol style="list-style-type: none"> <li>1. The laser is in STBY mode.</li> <li>2. The aiming light under menu is turned off.</li> <li>3. The fiber optic cable might be damaged or broken.</li> </ol>	<b>CAUSE:</b>	<ol style="list-style-type: none"> <li>1. Unit is not in ACTIVE MOVE.</li> <li>2. I am ACTIVE but my Activation Pedal is not working.</li> <li>3. Why is my Activation pedal vibrating?</li> </ol>
<b>SOLUTION:</b>	<ol style="list-style-type: none"> <li>1. Aiming light is only visible under ACTIVE mode. Select ACTIVE in the Guided Touch Interface.</li> <li>2. Select menu, then aiming light. With arrows UP/DOWN adjust the aiming light intensity.</li> </ol>	<b>SOLUTION:</b>	<ol style="list-style-type: none"> <li>1. Make sure you select ACTIVE and tap the Activation Pedal once to wake it up.</li> <li>2. Even under ACTIVE mode, you need to tap once at the pedal to wake it up.</li> <li>3. Gemini EVO laser has a Haptic Sense "HS" feature that enables you to sense small vibrations when the Activation Pedal is active. Under Menu and Haptic Sense "HS" you can select the intensity of the vibrations or completely turn it off.</li> </ol>

# SERVICE AND TROUBLESHOOTING

## TROUBLESHOOTING GUIDE

Why is my laser not connecting to Wi-Fi?	Photobiomodulation settings are not working and sounds strange
<p>CAUSE:</p> <ol style="list-style-type: none"><li>1. The unit has never been connected to a wireless network before.</li><li>2. The password my local network changed and I can't connect to the Wi-Fi.</li><li>3. I am not able to find my network while connecting to the Wi-Fi via the App.</li></ol>	<p>CAUSE:</p> <ol style="list-style-type: none"><li>1. After selecting treatment time, nothing happens.</li><li>2. After selecting a treatment time, the Activation Pedal is not activating the start of the procedure.</li><li>3. The unit sounds strange when using the 3 mm and 7 mm tips.</li></ol>
<p>SOLUTION:</p> <ol style="list-style-type: none"><li>1. Download the Gemini EVO laser mobile app and follow the instructions to make the Wi-Fi connection between your local network and your Gemini EVO device.</li><li>2. Under the Gemini EVO App, select "Add a Device" and follow the instructions. This will enable you to connect with your current network again with a new password.</li><li>3. Contact your network administrator to make sure your local network is broadcasting the Wi-Fi name. Sometimes for security purposes this feature is enabled, preventing outside users to see your local network.</li></ol>	<p>SOLUTION:</p> <ol style="list-style-type: none"><li>1. Once treatment time in seconds is selected, you need to select ACTIVE once more to enable the procedure.</li><li>2. If your Activation Pedal is not awake, tap once to turn it on. Then it will activate the procedure and enable multiple procedures consecutively.</li><li>3. Due to the low pulse modulation of power for the 3 mm and 7 mm tips, it is perfectly normal to hear the laser emit a deep pulse sound.</li></ol>

The display is dim or not illuminated?	General Questions
<p>CAUSE:</p> <ol style="list-style-type: none"><li>1. My unit display is barely visible but I can see some of the icons on the top part of the display.</li><li>2. Blank display: The unit lost power.</li><li>3. Only part of my display is visible.</li></ol>	<p>For a complete set of troubleshooting questions, videos and images on how your Gemini EVO laser should operate, please go to <a href="http://dashboard.geminievo.com">dashboard.geminievo.com</a> and select the troubleshooting tab. Additional information on setup, usability and settings can be found on <a href="http://dashboard.geminievo.com">dashboard.geminievo.com</a></p>
<p>SOLUTION:</p> <ol style="list-style-type: none"><li>1. It is likely your unit went into sleep mode. Simply touch anywhere on the Guided Touch Interface to wake up your Gemini Evo laser.</li><li>2. Power the unit on. The unit may have lost power due to low battery in which case the AC power supply should be plugged in.</li><li>3. The unit needs to be sent back to the manufacturer for repair. Contact your distributor representative for return instructions.</li></ol>	<p>For technical support and live troubleshooting with a technician, please contact our equipment support team at <a href="mailto:equipment.repair@ultradent.com">equipment.repair@ultradent.com</a> or 1.801.553.4574.</p>

## ERROR MESSAGES

-  **Software Update Error**  
The GEMINI EVO 810+980 Soft Tissue Laser is designed to perform periodic software updates. If during an update Internet connection is lost or unstable, the update may fail. The 'UF' error message shows on the display and user can restart laser unit to re-establish connectivity and resume the update.

---

-  **Overheating**  
The Gemini EVO 810+980 soft tissue laser is designed to perform surgical procedures at a specific temperature. High power and long procedures may cause the Laser Unit to heat up to the temperature threshold.  

Please wait a few minutes for the temperature to decrease before resuming normal operations.

---

-  **Activation Pedal Disconnected**  
The Gemini EVO 810+980 soft tissue laser is equipped with a long range Bluetooth chip.  

Please check the two AA batteries in the Activation Pedal and replace if needed. Press the Activation Pedal once to reactivate the connection with the Laser Unit. The Bluetooth icon on the Activation Pedal will turn blue, and the Bluetooth symbol will appear on the display when the laser is in Active mode and the Activation Pedal is successfully connected.

---

-  **Display Communication Error**  
There will be audible sound "Display Communication Error" should the glass electroluminescent display fail to turn on.  

Please plug the AC/DC power supply into the Laser Unit and restart the system by pressing the ON/OFF button. If the problem persists, contact the technical support for assistance.

---

-  **Calibration Error**  
The Gemini EVO laser is capable of sensing the internal laser light with a photodetector. If for any reason your Gemini EVO device goes out of calibration range, the CE ERROR message will be visible. At this time, we recommend contacting our technical support team as the unit might need to be sent for calibration.

# CYBERSECURITY

## OVERVIEW AND RECOMMENDATIONS

The Gemini EVO 810+980 soft tissue laser has been developed with Cybersecurity integrated throughout the total-product-lifecycle. Activities such as threat modeling, requirements documentation, penetration testing, and post market management planning have been executed for the device.

The Gemini EVO 810+980 soft tissue laser has been developed with Cybersecurity capabilities such as secure boot and code signing, using industry standard algorithms.

The Gemini EVO 810+980 soft tissue laser supports the ability to provide Cybersecurity Routine Updates and Patches remotely. The device provides notification on the Guided Touch Interface, mobile app and and web interface (Dashboard) when a new update is available. The user then has the ability to install the update directly into the device with any of these options. (Figure: 20.1)

A manufacturers statement on medical device security (MDS2) is available upon request for the Gemini EVO 810+980 soft tissue laser.

# ELECTROMAGNETIC ENVIRONMENT GUIDANCE

## ELECTROMAGNETIC COMPATIBILITY

**Notice:** The Gemini EVO 810+980 soft tissue laser complies with all requirements for electromagnetic compatibility according to IEC 60601-1-2: 2014.

**CAUTION:** Medical electrical equipment needs special precautions regarding electromagnetic compatibility (EMC) and needs to be installed and put into service according to the EMC information provided in the following tables. Portable and mobile Radio Frequency (RF) communications equipment can affect medical electrical equipment.

**WARNING:** The use of accessories, other than those specified, except those supplied or sold by Ultradent Products, Inc., as replacement parts for internal or external components, may result in increased EMISSIONS or decreased IMMUNITY of the Gemini EVO 810+980 soft tissue laser.

**Accessories:** Medical grade power supply - Maximum length 6ft (1.8 meters)

**Activation Pedal:** Wireless Bluetooth at 2.4GHz

Description: The Activation Pedal uses Bluetooth BLE 4.0 technology, which operates at a frequency of 2402 to 2480 MHz with TX power of +0dBm and RX sensitivity of -94.8dBm and uses GFSK modulation. The pedal is pre-configured by the manufacturer to only sync with the Gemini EVO Laser Unit that has a matching unique identifier. This prevents interference with other RF wireless technologies that may be present.

As a safety measure, any termination of the Bluetooth link between the Activation Pedal and the Laser Unit during use will result in the immediate termination of any laser emission. Reference the Service and Troubleshooting section of this manual should you encounter any connectivity issues between the Laser Unit and the Activation Pedal.

This device has passed wireless coexistence testing with common devices found in dental practices at a minimum separation distance of 30 cm.

## DEFINITIONS

**Emission (electromagnetic):** When electromagnetic energy is emitted by a source.

**Interference Immunity:** The ability of a device or system to work without errors even if there is electromagnetic interference.

**Immunity Level:** The maximum level of a certain electromagnetic interference that affects a particular device or system, where the device or system remains operative with a certain level of performance.

# ELECTROMAGNETIC ENVIRONMENT GUIDANCE

## ELECTROMAGNETIC EMISSION

The Gemini EVO 810+980 soft tissue laser is intended for operation in the electromagnetic environment specified below. The customer or user of the Gemini EVO laser should make sure that it is used in such an environment.

EMISSION TEST	COMPLIANCE	ELECTROMAGNETIC ENVIRONMENT – GUIDANCE
RF emissions according to CISPR 11	GROUP 1	The Gemini EVO laser uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions according to CISPR 11	CLASS B	The Gemini EVO laser is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions according to IEC 61000-3-2	CLASS A	
Voltage fluctuations/flicker emissions according to IEC 61000-3-3	COMPLIES	

## INTERFERENCE IMMUNITY

The Gemini EVO laser is intended for operation in the electromagnetic environment specified below. The customer or user of the Gemini EVO laser should make sure that it is used in such an environment.

INTERFERENCE IMMUNITY TEST	IEC 60601-1-2 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT GUIDANCE
Electrostatic discharge (ESD) according to IEC 61000-4-2	+/- 2, +/-4, +/- 6, +/-8 kV contact discharge +/- 2, +/-4, +/-8 kV, +/-15 kV air discharge	+/- 2, +/-4, +/- 6, +/-8 kV contact discharge +/- 2, +/-4, +/-8 kV, +/-15 kV air discharge	Floors should be wood, concrete, or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 50%.
Electrical fast transient/burst according to IEC 61000-4-4	100kHz repetition +/- 2kV for power supply lines 100kHz repetition +/- 2kV for power input/output lines	100kHz repetition +/- 2kV for power supply lines Not applicable	The quality of the line power supply should be that of a typical commercial or hospital environment.
Surge voltages according to IEC 61000-4-5	± 1 kV differential mode ± 2 kV common mode voltage	± 1 kV differential mode ± 2 kV common mode voltage	The quality of the line power supply should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and variations of the power supply according to IEC 61000-4-11	Voltage (VAC) 100, 240 over 0% residual with ½ cycles on a phase angle of 0° / 90° / 180° / 270°  Voltage (VAC) 100, 240 over 0% residual with 1 cycles on a phase angle of 0° / 90° / 180° / 270°  Voltage (VAC) 100, 240 over 70% residual with 25 cycles on a phase angle of 0° / 90° / 180° / 270°  Voltage (VAC) 100, 240 over 0% residual with 250 cycles on a phase angle of 0° / 90° / 180° / 270°	Voltage (VAC) 100, 240 over 0% residual with ½ cycles on a phase angle of 0° / 90° / 180° / 270°  Voltage (VAC) 100, 240 over 0% residual with 1 cycles on a phase angle of 0° / 90° / 180° / 270°  Voltage (VAC) 100, 240 over 70% residual with 25 cycles on a phase angle of 0° / 90° / 180° / 270°  Voltage (VAC) 100, 240 over 0% residual with 250 cycles on a phase angle of 0° / 90° / 180° / 270°	The quality of the line power supply should be that of a typical commercial or hospital environment.  If the user of the Gemini EVO laser requires it to continue functioning following interruptions of the power supply, it is recommended to have the Gemini EVO laser powered by an uninterruptible power supply or a battery.
Magnetic field of power frequencies (50/60 Hz) according to IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment

IMMUNITY TEST	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT GUIDANCE
Conducted RF IEC 61000-4-6  Radiated RF IEC 61000-4-3	3 Vrms 150 kHz to 80 MHz  3 V/m 80 MHz to 2.5 GHz	3 Vrms  3 V/m	Portable and mobile radio equipment must not be used within the recommended working clearance from the Gemini EVO laser and its cables, which is calculated based on the equation suitable for the relevant transmission frequency.  Recommended separation distance $d = [1.2] \sqrt{P}$ $d = [1.2] \sqrt{P}$ at 80 MHz to 800 MHz $d = [2.3] \sqrt{P}$ at 800 MHz to 2.5 GHz  where P is the nominal transmitter output in watts (W) specified by the transmitter manufacturer and d is the recommended working clearance in meters (m).  Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey <sup>2</sup> should be less than the compliance level <sup>3</sup> in each frequency range.  Interference is possible in the vicinity of equipment bearing the following graphic symbol.  

# ELECTROMAGNETIC ENVIRONMENT GUIDANCE

## NOTES

The higher frequency range applies at 80 MHz and 800 MHz.

Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast, cannot be predicted theoretically with accuracy. An investigation of the location is recommended to determine the electromagnetic environment resulting from stationary HF transmitters. If the measured field strength in the location in which the Gemini EVO laser is used exceeds the applicable RF compliance level above, the Gemini EVO laser should be observed to verify normal operation. If unusual performance characteristics are observed, it may be necessary to take additional measures such as reorientation or repositioning of the Gemini EVO laser.

Over the frequency range 150kHz to 80MHz, field strengths should be less than 3V/m.

## WORKING CLEARANCES

The Gemini EVO laser is intended for operation in an electromagnetic environment where radiated HF interference is checked. The customer or the user of the Gemini EVO laser can help prevent electromagnetic interference by duly observing the minimum distances between portable and/or mobile RF communication devices (transmitters) and the Gemini EVO laser. These values may vary according to the output power of the relevant communication device as specified below.

RATED MAXIMUM OUTPUT POWER OF TRANSMITTER [W]	WORKING CLEARANCE ACCORDING TO TRANSMISSION FREQUENCY [M]		
	150 KHZ TO 80 MHZ	80 MHZ TO 800 MHZ	800 MHZ TO 2.5 GHZ
	$d = [1.2] \sqrt{P}$	$d = [1.2] \sqrt{P}$	$d = [2.3] \sqrt{P}$
0.01	0.12	0.12	.023
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters whose maximum nominal output is not specified in the above table, the recommended working clearance in meters (m) can be determined using the equation in the corresponding column, where P is the maximum nominal output of the transmitter in watts (W) specified by the transmitter manufacturer.

**Remark 1:** The higher frequency range applies at 80 MHz and 800 MHz.

**Remark 2:** These guidelines may not be applicable in all cases. The propagation of electromagnetic waves is influenced by their absorption and reflection by buildings, objects, and persons.

## EMERGENCY TERMINATION OF LASER EMISSIONS

The Gemini EVO 810+980 soft tissue laser has been designed with several methods to terminate emission of laser energy in emergency situations. These methods include a power button (ON/OFF) (Figure: 21.1) and the emergency (STOP) (Figure: 21.2) button located at the front of the laser unit.

## WARRANTY

Your Gemini EVO laser comes with a 2-year factory warranty. Extended warranty can be purchased and warranty certificates can be downloaded by accessing your Dashboard page at [dashboard.geminievo.com](http://dashboard.geminievo.com) or under the warranty icon in the iOS and Android Apps.

Seller warrants the Products to be free from defects in materials and workmanship for a period of twenty-four months from the date of shipment, except for consumables. If within such period any Products shall be proved to Seller's satisfaction to be defective, it shall be (i) repaired using new or refurbished parts, or (ii) replaced with a new or refurbished product, at Seller's sole discretion. Such repair or replacement shall be Seller's sole obligation and Buyer's exclusive remedy under this Warranty and shall be conditioned, at Seller's option, upon return of such Products to Seller, f.o.b. its factory. This Warranty only covers Product issues caused by defects in material or workmanship during ordinary consumer use; it does not cover Product issues caused by any other reason, including but not limited to acts of God, modifications of or to any part of the Product, improper testing, assembly, mishandling, misuse, neglect, adjustments, alterations to the products, improper operation contrary to current instructions relating to installation, maintenance or operation, or contrary to industry standards relating to acceptable input power.

THIS WARRANTY IS EXCLUSIVE AND IN LIEU OF ALL OTHER REPRESENTATIONS AND WARRANTIES, EXPRESS OR IMPLIED; AND SELLER EXPRESSLY DISCLAIMS AND EXCLUDES ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE. SELLER SHALL HAVE NO OBLIGATION OR LIABILITY TO REFUND ANY PORTION OF THE PURCHASE PRICE AND SHALL NOT BE LIABLE FOR ANY SPECIAL, EXEMPLARY, INCIDENTAL, CONSEQUENTIAL OR PUNITIVE DAMAGES OR LOST PROFITS, OR DAMAGE TO PERSON OR INJURY IN CONNECTION WITH THE PURCHASE OR USE OF THE INSTRUMENT.

SYMBOLS	DESCRIPTION	SYMBOLS	DESCRIPTION
	<b>MANUFACTURER</b> INDICATES WHICH COMPANY MANUFACTURES		<b>SHIP VERTICAL, WITH ARROWS POINTED UPWARD</b>
	<b>DATE OF MANUFACTURE</b> INDICATES THE DATE AND YEAR OF MANUFACTURE		<b>FRAGILE - HANDLE WITH CARE</b>
	<b>CATALOG PART NUMBER</b> INDICATES THE MANUFACTURER PART NUMBER.		<b>DO NOT USE IF PACKAGE IS DAMAGED</b>
	<b>SERIAL NUMBER</b> INDICATES SERIAL NUMBER FOR THE PRODUCT PART		<b>RECOMMENDED STORAGE TEMPERATURE</b>
	<b>PRESCRIPTION STATEMENT</b> FEDERAL LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A DENTIST OR PHYSICIAN OR OTHER LICENSED MEDICAL PRACTITIONER		<b>ATMOSPHERIC PRESSURE LIMITATION</b>
	<b>LASER APERTURE</b> INDICATES THE SYSTEM CONTAINS A LASER		<b>RELATIVE HUMIDITY RANGE</b>
	<b>WARNING</b> INDICATES POSSIBLE EXPOSURE TO BOTH RED AND INFRARED LASER RADIATION		<b>KEEP AWAY FROM HEAT/SUNLIGHT</b>
	<b>OPTICAL FIBER APPLICATOR</b> INDICATES WHERE LASER ENERGY COMES OUT		<b>KEEP DRY</b>
UNIT LABEL	PEDAL LABEL	SYMBOLS	DESCRIPTION
			<b>TYPE B APPLIED PART</b> THE APPLIED PART IS NOT CONDUCTIVE TO THE PATIENT
			<b>MUST REFER TO USER'S MANUAL</b>
			<b>LASER STOP</b> EMERGENCY SWITCH TO STOP LASER OUTPUT POWER
			<b>NON-IONIZING RADIATION</b>
	<b>HANDPIECE</b> OPTICAL FIBER		<b>WESTERN EUROPE ELECTRONIC AND ELECTRICAL DISPOSAL STANDARD (WEEE)</b>

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Nominal Ocular Hazard Distance (NOHD) is the distance from the source of laser emission to the point where it no longer exceeds its Maximum Permissible Exposure (MPE – highest level of laser radiation to which a person may be exposed without hazardous effects or adverse biological changes in the eyes or skin). The Nominal Hazard Zone (NHZ) is the space within which the level of direct, reflected, or scattered radiation during normal operation exceeds the appropriate MPEs.

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Manufactured by Azena Medical, LLC. 3021 Citrus Circle, Suite 180, Walnut Creek, CA 94598 USA.
Made in the USA from U.S. and globally sourced materials. Distributed by Ultradent Products Inc.,
**MedEnvoy Global B.V.**