TABLE OF CONTENT

bg	Инструкции за употреба	4
cs	Návod k obsluze	13
da	Betjeningsvejledning	20
de	Gebrauchsanweisung	27
el	Οδηγίες Χρήσης	35
en	Instructions for use	43
es	Instrucciones de uso	50
fi	Käyttöohje	58
fr	Notice d'utilisation	65
hr	Upute za uporabu	73
hu	Használati útmutató	80
it	Istruzioni per l'uso	87
lt	Naudojimo instrukcija	94
lv	Lietošanas pamācība	101
nl	Gebruiksaanwijzing	108
no	Bruksanvisning	116
pl	Instrukcja obsługi	123
pt	Instruções de utilização	131
ro	Instrucţiuni de utilizare	139
sk	Návod na použitie	146
sl	Navodila za uporabo	153
sr	Uputstvo za upotrebu	160
sv	Bruksanvisning	167
tr	Kullanım talimatları	174

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Zimmerman EE, Crawford P. Cutaneous cryosurgery. Am Fam Physician. 2012 Dec 15;86(12):1118-24.
Clebak KT, Mendez-Miller M, Croad J. Cutaneous Cryosurgery for Common Skin Conditions. Am Fam Physician. 2020 Apr 1;101(7):399-406.
Pasquali, Paola. (2015). Cryosurgery: A Practical Manual. 10.1007/978-3-662-43939-5.
Usatine, R.P., Stulberg, D.L., & Colver, G.B. (2005). Cutaneous Cryosurgery: Principles and Clinical Practice (4th ed.). CRC Press.
<https://doi.org/10.1201/b17660>

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EN – CRYOPRO® INSTRUCTIONS FOR USE

TABLE OF CONTENT

Intended use	Page 43
Adverse effects.....	Page 43
Intended Target Groups.....	Page 43
Indications.....	Page 43
Contraindications.....	Page 43
Warnings.....	Page 43
Filling the unit.....	Page 44
Accessories.....	Page 44
Operating Principles.....	Page 44
Treatment/Clinical application	Page 44
Decontamination	Page 44
Storing	Page 45
Maintenance	Page 45
Disposal	Page 45
Performance Characteristics	Page 45
Explanation of symbols.....	Page 45
Accessories.....	Page 46
Recommended treatment times.....	Page 47
Cleaning and sterilization of CryoPro® contact probes	Page 48

The CryoPro is an active non-implantable medical device for cryosurgical treatment of skin lesions by tissue freezing. The CryoPro device is intended to be filled with liquid nitrogen (LN). On activation, the LN is expelled through a replaceable filter and an immersion tube to the tip of the device where various spray apertures or contact probes can be attached.

INTENDED USE

The CryoPro® is a handheld device for cutaneous cryosurgery that uses localized application of liquid nitrogen to treat skin lesions by tissue freezing. Treatment can be performed using open spray or contact probe techniques by healthcare professionals such as dermatologists, family physicians, and nurses trained in cryosurgery in a normal well-lit and ventilated clinical setting.

The use of cryosurgery is well tolerated, fast, easy to apply, and without systemic side effects.

ADVERSE EFFECTS

Temporary adverse effects;

- Pain during and/or shortly after treatment, bulla formation and local edema are the major temporary adverse effects.

Long term complications;

- In some cases, the patient may develop lesional hypopigmentation and/or peripheral hyperpigmentation. Risk of permanent hypopigmentation is higher in dark-skinned persons.

INTENDED TARGET GROUPS

The device can be used for all groups of patients regardless of age, gender and race. However, infants as well as dark skinned people not willing to accept pigment changes, are excluded.

INDICATIONS

The CryoPro® is indicated for the following benign skin lesions:

Acne, dermatofibroma, hemangioma, hypertrophic scar, keloid, seborrheic keratosis, myxoid cyst, pyogenic granuloma, skin tag, solar lentigo, and viral warts (anogenital, common, flat, mosaic, plantar).

Several premalignant/malignant lesions may also be treated using cryotherapy:

Actinic keratosis, basal cell carcinoma (BCC), Bowen's disease, squamous cell carcinoma (SCC).

The CryoPro® should NOT be used outside of the above-mentioned indications.

CONTRAINDICATIONS

- BCC/SCC lesions with high risk features
- Lesions where pathology is needed
- Previously experienced sensitivity to cryosurgery
- Immunosuppression
- Impaired vascular supply
- Anticoagulant use
- Treatment on areas above prominent bones or superficial nerves
- Tumors with indefinite margins
- Not willing to accept pigment changes (higher risk in dark-skinned people)
- Infants
- For patients with psoriatic arthritis, cryosurgical treatment may in rare cases be a trigger for flaring.

For additional information about contra-indications, please see relevant literature, e.g. Clebak KT et al. Cutaneous Cryosurgery for Common Skin Conditions. Am Fam Physician. 2020 Apr 1;101(7):399-406.

WARNINGS

- Do never use the CryoPro® unit without either a spray aperture or contact probe in place. The excessive amount of liquid put out if the unit is used without an attached spray or probe may cause severe damage to the patient. A Flow Restrictor (accessory) can be affixed to prevent nitrogen flow if a spray aperture or contact probe is not in place.
- After filling, the lid must be mounted and tightened securely.
- The CryoPro® is supposed to be held in a vertical position when used. Liquid nitrogen may escape through the pressure relief valve if the unit is held in a position allowing the liquid inside the bottle to be in contact with the lid. The burst of outflowing nitrogen may hit the hand of the user and/or the patient.
- Do not use the unit if the pressure relief valve is blocked or have been tampered with or the unit is otherwise damaged. This may cause excessive pressure and violent release of nitrogen, exposing both patient, user and other persons to danger.
- The unit may get cold (visible ice on lid and spraytube) after prolonged use. Leave the unit to thaw before resuming treatment.
- Be aware of signs of insulation failure – the bottle will freeze over. DO NOT USE if the bottle is covered with ice.
- Be careful when opening a pressurized unit. Open the lid slowly. The unit may still be pressurized even if no residual liquid is left.
- DO NOT spray on open wounds, into body orifices, or into the eyes.
- Only use liquid nitrogen in a well ventilated room. Evaporated nitrogen may cause asphyxiation.
- Do not use contact probes with signs of a compromised surface. Scratches or loosening of the surface coating may bring the patient into contact with underlying materials or debris.
- Do not use the CryoPro® unit if there are signs of severe transport damage or unintentional opening of the packaging before first use.
- Gloves, goggles, and suitable protective clothing should be worn during the cleaning process.
- Failure to remove contaminants from the surface of the contact probe could compromise the correct functioning, its safe use, and the subsequent sterilization process.
- The CryoPro® is a medical device intended for professional use only. When not in use it must be kept out of reach of children and unauthorized personnel.
- When using a contact probe, the sili-cone vent tube, which is attached to the base of the contact probe, must be directed away from the patient or user until hardened. This will happen approximately 5 seconds after activating the liquid nitrogen flow and the tube will remain in a fixed position (frozen) during the remainder of the procedure.

NB: Any serious incident that occurs in relation to the device must be reported to the manufacturer and the relevant competent authority of the EU member state in which the user and/or patient is established. For users outside the EU: local distributor and relevant competent authorities must be informed.

FILLING THE UNIT

The CryoPro® unit is to be filled with *liquid nitrogen only*. Fill the unit by slowly pouring liquid nitrogen into the bottle or by using a standard low pressure withdrawal device from a storage tank (dewar). Depending on the number of lesions to be treated and the duration of the individual treatments the unit should be filled to at least one third of its maximum capacity. After filling a warm unit, the liquid will boil for approx. 30 secs. Hereafter, the lid can be mounted and tightened firmly. A weak but audible hissing sound from the pressure relief valve indicates that normal working pressure has been reached inside the unit. This is normal operation and the pressure relief valve should never be blocked or tightened. Various factors influence the time it takes from the lid is closed until the device is ready for use. Please see recommendations in the section "Performance characteristics".

Caution should be taken when refilling a CryoPro® containing residual fluid as the unit is pressurized at 0.8 bar. Slowly unscrewing the lid will allow the unit to depressurize before the top comes off. This is accompanied by an audible hissing sound until it is fully depressurized. Evaporized nitrogen may escape and hit the hand.

NB: When handling liquid nitrogen (LN₂) it is recommended to take appropriate precautions as prescribed in the Safety Data Sheet (SDS) from the LN₂-supplier. As a minimum, it is recommended to use protective gloves and safety goggles or similar.

ACCESSORIES

Before use, attach a spray aperture or contact probe, as the CryoPro cannot perform its intended use without an attached spray aperture or contact probe.

The CryoPro® is supplied with a set of standard open spray apertures. Optional accessories comprise a variety of spray tips, contact probes, rubber/plastic covers (cones), a replaceable filter and other attachments. A table listing the accessories are found in the section "Accessories".

NB: The contact probes are semi-critical items that come into contact with non-intact skin. It is important to inspect them as they must be discarded if they have any signs of a compromised surface. Scratches or defects in the gold coating may bring the patient into contact with underlying materials or debris.

NB: The CryoPro® is only intended for use with Cortex manufactured spray tips, contact probes, and accessories.

OPERATING PRINCIPLE

When the device is pressurized and the trigger handle is activated, the flow of LN₂ will start. The nitrogen flow will stop when the trigger handle is released. Adjustment of the LN₂ flow rate is not possible using the trigger handle. The nitrogen flow rate is adjusted by selecting the appropriate spray aperture. A continuous flow of LN₂ without pressing the handle can be achieved by pushing the handle forward and to the side. For left-handed users, the trigger handle can be placed on the left side of the main valve.

A threaded replaceable filter is attached to the inlet of the immersion tube, to prevent clogging of tube, valve and spray tip's. The filter is an accessory and can be replaced if needed. Remove and tighten the filter by hand only. Do not use tools. Make sure the filter is properly seated before mounting the lid on the bottle.

NB: To avoid contamination and cold burn, it is recommended to use protective gloves and safety goggles or similar.

TREATMENT/CLINICAL APPLICATION

All probes and spray apertures are mounted into the permanently affixed knurled nut with fingertight firmness. Do not apply excessive force.

Anesthetic cream may be applied to improve tolerability of the procedure and reduce posttreatment pain.

NB: When left on the table between uses, there will be an audible, constant, or intermittent hissing sound from the pressure relief valve. This is normal and should not be prevented.

Spray freezing (REF OS A, B, C, D, SS, BS, SOFT, LL)

The spray apertures include: A (1.00 mm opening), B (0.75 mm opening), C (0.55 mm opening), D (0.45 mm opening) and conical spray apertures: A (1.00 mm opening), B (0.75 mm opening), C (0.55 mm opening), D (0.45 mm opening), as well as Bent Spray apertures are available: BS BL (75mm long 0.75 mm opening) and (BS CL, (75mm long 0.55 mm opening) and four Straight Spray apertures: (SS CL, (75mm long 0.55 mm opening), SS AS (25mm long 1.00 mm opening), SS BS (25mm long 0.75 mm opening), SS CS (25mm long 0.55 mm opening), both all for areas difficult to reach like behind the ears, lip and outer nose (to avoid nitrogen gas into the airways).

Select a spray aperture providing an LN₂ output appropriate for the lesion to be treated. In-depth freezing is best obtained when the spray aperture is close to the lesion (5 - 10 mm distance) as the liquid content of the spray is higher close to the aperture. When possible, raised warts should be approached tangentially from at least two sides. This allows the ice ball to move down through the wart while sparing the surrounding tissue. Spraying directly at the lesion (perpendicular to the surface) will cause more lateral spread with less penetration. To avoid this apply the spray intermittently or use a smaller spray aperture.

The Soft Spray aperture is for freezing larger areas at greater distance (up to approx. 5 cm) with higher vapor content. Move the aperture slowly back and forth thus "painting" the surface with vaporized nitrogen.

The Luer Lock adapter facilitates the use of a hypodermic needle as a spray aperture and the diameter of the needle determines the LN₂ flow rate. The flow rate must be adjusted according to the lesion to be treated.

NB: The use of hyperdermic needles is not intended for intralesional treatment.

If the spray is clogged or liquid nitrogen flows slowly, it might be due to an impurity in the liquid nitrogen. Take off the lid and let it thaw completely. Inspect the filter and replace if necessary. Only use original Cortex replacement filter. To clean the spray aperture, use the cleaning adapter or an air compressor.

Contact freezing (REF CP SP, CP 1mm-30mm)

Closed contact probes provide in-depth freezing with minimal lateral spread.

Select a probe size fitting the size of the lesion. Make the lesion wet with a drop of water or contact gel prior to freezing to facilitate the thermal transfer between skin and probe tip. Apply the contact probe when still warm applying light pressure and pull the trigger handle. Keep freezing until a 1-5 mm halo is seen around the probe tip. Then stop the freezing and allow a few seconds for the probe tip to be released from the skin.

Freeze time varies according to the aperture or probe used, the size and type of the lesion and the freezing distance (spray only). The user is encouraged to familiarize him/her-self with the available literature on cryosurgery using liquid nitrogen. A table showing the recommended treatment times for the indicated skin lesions is seen in the section "Recommended treatment times".

After treatment, the treated area must be protected (normal wound care).

DECONTAMINATION

The main CryoPro® unit is not intended to be autoclaved, but may be cleaned using alcohol (isopropyl) or ethanol). Normally, spray apertures need no cleaning as they are not in contact with the patient.

The CryoPro® contact probes are semi-critical reusable devices that may get in contact with non-intact skin during the intended use. For this reason, a manual cleaning process together with a moist heat sterilization process have been validated as being able to prepare the device for reuse. See section "Cleaning and sterilization of CryoPro® contact probes".

STORING

Storing liquid nitrogen

To ensure a clean supply of liquid nitrogen the storage tank (dewar) must be completely emptied prior to refilling 3 - 4 times per year. If necessary, water from icecrystals remaining inside the tank can be dried out by creating air circulation (e.g. using a vacuum cleaner for approx. 30 minutes).

Storing the CryoPro® device

At the end of the workday it is highly recommended to empty the unit by removing the lid and discarding any excess liquid nitrogen. Do NOT pour it back into the dewar to avoid contamination. Once the bottle is empty it is recommended to store the bottle upside down for a minimum of 30 minutes before applying the lid. In order to avoid accumulation of condensation the CryoPro® should be stored with the lid on.

The empty device and accessories must be stored at regular room temperature and humidity. There are no special storage requirements and no restrictions on shelf-life. Contact probes do not need to be stored in sterile condition but must be sterilized before use.

The lifetime of the device and the accessories is 5 years. The lifetime of the contact probes is 5 years or 100 reprocessing cycles.

MAINTENANCE

No preventive maintenance is needed.

DISPOSAL

CryoPro® and metal accessories can be disposed of as metal scrap or returned to Cortex Technology for disposal. It is recommended to decontaminate contact probes before disposal or dispose as contaminated materials.

PERFORMANCE CHARACTERISTICS

Device	LN ₂ capacity	Static holding time
CryoPro® MAXI	500 ml	20-24 hours
CryoPro® MINI	350 ml	12-14 hours

The flow rate and the velocity of the LN₂ and the maximum treatment time depends on the spray tip (A-D) mounted:

Spray tip	A	B	C	D
Hole diameter (mm)	1.00	0.75	0.55	0.45
Average flow within the first minute (mL/min)	54	28	15	9
Average velocity (m/s)	0.75	0.88	0.92	0.98
Estimated total treatment time in 1 min sessions (MINI/MAXI) (minutes)	6/9	12/17	22/32	36/52

Approximate time from closing the lid after filling the CryoPro® until it is ready for use			
CryoPro® thermal state when filled	Room temperature (22±1 °C)		Cooled (refilled)
Fill level	Maximum volume	1/3 of max. volume	1/3 of max. volume
CryoPro® Maxi (500 mL)	30 seconds	1 minute and 30 seconds	15 minutes
CryoPro® Mini (350 mL)	25 seconds	1 minute	10 minutes

EXPLANATION OF SYMBOLS

Graphic	Description
	Indicates the medical device manufacturer.
	Indicates the date when the medical device was manufactured.
	Indicates the manufacturer's catalogue number so that the medical device can be identified.
	Indicates the manufacturer's serial number so that a specific medical device can be identified.
	Indicates the need for user to consult the instructions for use.
	Indicates the item is a medical device.
	Indicates a carrier that contains unique device identifier information.
	CE-marking with Notified Body reference no. xxxx. Signifies European technical conformity.
	Caution: Federal (US) law restricts this device to sale by or on the order of a physician.

ACCESSORIES

All accessories are compatible with both CryoPro® Maxi and Mini.

Spray tips	Classification	Ref. code:	Purpose:
Spray Tip "A", 1,00mm	Ila	OS A	Open spray aperture (1.00 mm)
Spray Tip "B", 0.75 mm	Ila	OS B	Open spray aperture (0.75 mm)
Spray Tip "C", 0.55 mm	Ila	OS C	Open spray aperture (0.55 mm)
Spray Tip "D", 0.45 mm	Ila	OS D	Open spray aperture (0.45 mm)
Conical Spray Tip "A", 1.00 mm	Ila	COS A	Open spray aperture (1.00 mm)
Conical Spray Tip "B", 0.75 mm	Ila	COS B	Open spray aperture (0.75 mm)
Conical Spray Tip "C", 0.55 mm	Ila	COS C	Open spray aperture (0.55 mm)
Conical Spray Tip "D", 0.45 mm	Ila	COS D	Open spray aperture (0.45 mm)
Bent Spray extension "C", 0.55 mm, 75 mm long	Ila	OS BS CL	Spray tip used for accessing hard to reach ana-tomical structures
Bent Spray extension "B", 0.75 mm, 75 mm long	Ila	OS BS BL	Spray tip used for accessing hard to reach ana-tomical structures
Straight Spray extension "C", 0.55 mm 75 mm long	Ila	OS SS CL	Spray tip used for accessing hard to reach ana-tomical structures
Straight Spray extension "A", 1.00 mm , 25 mm short	Ila	OS SS AS	Spray tip used for accessing hard to reach ana-tomical structures
Straight Spray extension "B", 0.75 mm, 25 mm short	Ila	OS SS BS	Spray tip used for accessing hard to reach ana-tomical structures
Straight Spray extension "C", 0.55 mm ,25 mm short	Ila	OS SS CS	Spray tip used for accessing hard to reach ana-tomical structures
Acne Soft Spray Tip	Ila	OS SOFT	Spray tip used for treating acne

Contact probes:	Classification	Ref. code:	Purpose:
Sharp Pointed Contact probe	Ila	CP SP	Contact probe
1 mm Contact probe	Ila	CP 1 mm	Contact probe
2 mm Contact probe	Ila	CP 2 mm	Contact probe
3mm Contact probe	Ila	CP 3 mm	Contact probe
4 mm Contact probe	Ila	CP 4 mm	Contact probe
5 mm Contact probe	Ila	CP 5 mm	Contact probe
6 mm Contact probe	Ila	CP 6 mm	Contact probe
8 mm Contact probe	Ila	CP 8 mm	Contact probe
10 mm Contact probe	Ila	CP 10 mm	Contact probe
15 mm Contact probe	Ila	CP 15 mm	Contact probe
20 mm Contact probe	Ila	CP 20 mm	Contact probe
30 mm Contact probe	Ila	CP 30 mm	Contact probe

Covers	Classification	Ref. code:	Purpose:
Rubber Cones	I	RUB	Used as a mechanical barrier to confine the open spray to affect only skin intended to be treated
Lexan Disc	I	LEXAN	Used as a mechanical barrier to confine the open spray to affect only skin intended to be treated

Attachments:	Classification	Ref. code:	Purpose:
Open Flow restrictor	N.A.	OFR	Blocks liquid nitrogen flow when no nozzle or contact probe is attached
Luer Lock adapter	Ila	OS LL	Allows hypodermic needle compliant with 'ISO 80369-7 female luer lock interface' to be attached
Spray Tip Cleaning adapter	N.A.	TCA	Used for flushing clogged spray apertures (re-moves particles which are blocking the flow)
Replaceable Nitrogen Filter	Ila	RNF	Replaceable filter for the inlet end of the immersion tube

For a visualization of the accessories: See enclosed "Z11006.xx Treatment accessories and handling devices" on our website www.cortex.dk.

RECOMMENDED TREATMENT TIMES

The table below shows a list of recommendations for treating various skin lesions using liquid nitrogen. The technique, freeze time in seconds, number of freeze-thaw-cycles (FTC), margin around the treated lesion, and number of treatment sessions are shown for indicated benign and premalignant/malignant skin lesions. The recommended cryosurgical technique can be either open spray (OS), contact probe (CP)

Benign lesions	Technique	Freeze Time (s)	FTC	Margin (mm)	Treatment sessions
Acne	OS/CP/SS	5 to 15	1	1	1
Dermatofibroma	OS/CP	20 to 60	1	2	2
Hemangioma	OS/CP	10	1	<1	1
Hypertrophic scar	OS/CP	20	1	2	1
Keloid	OS/CP	30	1-3	2	3
Seborrheic Keratosis	OS	10 to 15	1-3	<1	1 to 3
Myxoid Cyst	OS/CP	5	1	1 to 2	1
Pyogenic Granuloma	OS	15	1	<1	1
Skin Tags	OS/CP	5	1	1 to 2	1
Solar Lentigo	OS	5	1	<1	1
Viral Warts					
Anogenital (condyloma acuminata)	OS	10	1	2	3 to 4
Common (verruca vulgaris)	OS	10	1	1	2
Flat (verruca plana)	OS	10	1	1	2
Mosaic	OS	30	1	3 to 5	4
Plantar (verruca plantaris)	OS	20	1	2	3 to 4
Premalignant/ malignant lesions:	Technique	Freeze Time (s)	FTC	Margin (mm)	Treatment sessions
Actinic Keratosis	OS	5 to 20	1	1 to 2	1
Basal Cell Carcinoma (BCC)	OS	60 to 90	1 to 3	5	1
Bowen's disease	OS	15 to 30	1 to 2	3	1
Squamous Cell Carcino-ma (SCC)	OS	60 to 90	1 to 3	5	1

OS: Open Spray

CP: Contact Probe

SS: Soft Spray

The table above is based on recommendations in the literature:

Zimmerman EE, Crawford P. Cutaneous cryosurgery. Am Fam Physician. 2012 Dec 15;86(12):1118-24.

Clebak KT, Mendez-Miller M, Croad J. Cutaneous Cryosurgery for Common Skin Conditions. Am Fam Physician. 2020 Apr 1;101(7):399-406.

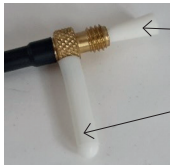
Pasquali, Paola. (2015). Cryosurgery: A Practical Manual. 10.1007/978-3-662-43939-5.

Usatine, R.P., Stulberg, D.L., & Colver, G.B. (2005). Cutaneous Cryosurgery: Principles and Clinical Practice (4th ed.). CRC Press.

<https://doi.org/10.1201/b17660>

CLEANING AND STERILIZATION OF CRYOPRO® CONTACT PROBES

According to ISO 17664-1:2021 Table B1 – Processing instructions (reusable medical devices).

Manufacturer: Cortex Technology Aps		Method:	Symbol:
Device(s): CryoPro® contact probes (REF CP SP, CP 1mm-30mm)			
WARNINGS	The devices are delivered in a non-sterile state and must be cleaned and sterilized before use.		
	Gloves, goggles, and suitable protective clothing should be worn during the cleaning process.		
	Failure to remove contaminants from the surface of the contact probe could com-promise the correct functioning, its safe use, and the subsequent sterilization process.		
	Only sterilization using hot steam is recommended. Other sterilization methods based on e.g. EtO gas, liquid chemicals, dry heat, VHP, or ionizing radiation have not been validated to be effective and may cause physical harm to the product.		
Limitations on processing:	To avoid burns, allow the devices too cool after sterilization (minimum 30 minutes).		
	Re-processing has minimal effect on the contact probes. The end of life is determined by physical wear and tear. Do not use the contact probes if the gold coated applied part is discolored or has signs of physical damage (e.g., scratches). The threaded brass part of the probe will turn brown/black as it tarnishes. This is a normal process that does not affect the functionality of the product and is of purely cosmetic concern.		
	The contact probes have successfully been validated to withstand up to 100 cleaning and sterilization cycles. Replace the probe after 100 reprocessing cycles or in case of discoloration or physical damage to the gold coating on the part that may get in contact with the patient.		
INSTRUCTIONS			
Initial treatment at the point of use	To avoid contamination, the contact probe should be kept in the sterilization pouch during transport from the autoclave to the point of use. When opening the sterilization pouch, be careful not to contaminate the probe before use. The protective silicone cap must be replaced by the silicone vent tube and the silicone stud must be removed before use.		
	At the point of use it must be inspected if the contact probe surface is intact. The contact probe must be discarded in case of discoloration or surface defects in the gold coating.		
	After use, detach the contact probe from the CryoPro®. Remove any gross soiling using a clean lint-free wipe.		
	Before re-use, the contact probe must be cleaned and sterilized according to the validated procedure described below.		
Preparation before cleaning	Replace the silicone vent tube on the contact probe with the protective silicone cap (see image below). Also, attach the silicone stud in the brass thread, as any liquid ingress may block the nitrogen flow during use when it freezes.		
			
Cleaning: Automated	The contact probes are not suitable for ultrasonic or automated cleaning processes due to the risk of liquid ingress that may block the nitrogen flow. Therefore, only a manual cleaning method is provided.		

INSTRUCTIONS													
Cleaning: Manual	Cleaning agent: Prolystica 2X Concentrate Neutral Detergent from the manufacturer Steris. Composition of the detergent is shown below:												
	<table><tr><th>Name</th><th>Product identifier</th></tr><tr><td>Ethoxylated Alcohol Surfactant</td><td>(CAS No) 68439-46-3</td></tr><tr><td>Poly(oxy-1, 2-ethaediyl), .alpha.-phenyl-.omega.-hydrixy-</td><td>(CAS No) 9004-78-8</td></tr><tr><td>1-Octanamine, N, N-dimethyl-, N-oxide</td><td>(CAS No) 2605-78-9 (REACH No) 01-2119409076-45-0002</td></tr><tr><td>Methyl-oxirane polymer with oxirane</td><td>(CAS No) 9003-11-6</td></tr></table>			Name	Product identifier	Ethoxylated Alcohol Surfactant	(CAS No) 68439-46-3	Poly(oxy-1, 2-ethaediyl), .alpha.-phenyl-.omega.-hydrixy-	(CAS No) 9004-78-8	1-Octanamine, N, N-dimethyl-, N-oxide	(CAS No) 2605-78-9 (REACH No) 01-2119409076-45-0002	Methyl-oxirane polymer with oxirane	(CAS No) 9003-11-6
	Name	Product identifier											
	Ethoxylated Alcohol Surfactant	(CAS No) 68439-46-3											
	Poly(oxy-1, 2-ethaediyl), .alpha.-phenyl-.omega.-hydrixy-	(CAS No) 9004-78-8											
	1-Octanamine, N, N-dimethyl-, N-oxide	(CAS No) 2605-78-9 (REACH No) 01-2119409076-45-0002											
	Methyl-oxirane polymer with oxirane	(CAS No) 9003-11-6											
Cleaning instructions: 1. Dilute the detergent in tap water (1 ml detergent per 1 liter of water, as worst-case dosage). 2. Once the protective silicone caps are mounted, immerse the entire contact probe in the solution and leave for 5 minutes. 3. Carefully remove the contact probe from the solution and brush the contact probe for minimum 1 minute to remove visible debris. Use a soft nylon brush that does not scratch the surface. 4. Rinse the contact probe under running water and brush until visibly clean for at least 1 minute. Tap water can be used. The cleaning procedure must be performed prior to sterilization.													
Disinfection	The CryoPro contact probes are not intended to be disinfected. A sterilization instruction is provided below.												
Drying	Dry the contact probe with a clean lint-free wipe and prepare for sterilization.												
Maintenance, inspection, and testing	1. Inspect the contact probe for any visible debris. 2. Inspect the contact probe for any signs of a compromised surface. The contact probe must be discarded in such case, as scratches or loosening of the surface coating may bring the patient into contact with underlaying materials or debris that are not biocompatible.												
Packaging	Place the cleaned contact probe in a single pouch suitable for steam sterilization (in accordance with EN ISO 11607-1:2020).												
Sterilization	Make sure that the silicone vent tube is not attached, and that the protective silicone cap and stud are in place during the sterilization process to prevent steam/water from entering the lumen inside the device. In case water is trapped inside the probe it will freeze during use and block the nitrogen flow. Processing parameters that have been validated to effectively reach a sterility as-surance level (SAL) of 1x10-6 are summarized below:												
	<table><tr><th>Autoclave cycle type</th><th>Temperature</th><th>Exposure time (Min.)</th><th>Drying time (Min.)</th></tr><tr><td>Gravity displacement</td><td>121 °C (250 °F)</td><td>15 minutes</td><td>20 minutes</td></tr></table>			Autoclave cycle type	Temperature	Exposure time (Min.)	Drying time (Min.)	Gravity displacement	121 °C (250 °F)	15 minutes	20 minutes		
	Autoclave cycle type	Temperature	Exposure time (Min.)	Drying time (Min.)									
Gravity displacement	121 °C (250 °F)	15 minutes	20 minutes										
Allow the sterilized devices to properly cool to room temperature in a dedicated area with minimum traffic and air flow for minimum 30 minutes. Do not open the sterilization pouch until the device has been transported to the point of use. The silicone vent tube should not be attached until the device has been transported to the point of use.													
Storage	CryoPro contact probes must be stored at regular room temperature and humidity. Apart from that, there are no special storage requirements and no restrictions on shelf-life. The contact probes are intended to be sterilized before use and can be stored according to the specifications of the sterilization pouch manufacturer.												
Additional information													
Manufacturer contact	Cortex Technology, Niels Jernes Vej 6B, 9220 Aalborg, Denmark Phone: +45 98574100, Email: cortex@cortex.dk												

The instructions provided above have been validated by the manufacturer of the medical device as being capable of preparing the contact probes for reuse. It remains the responsibility of the processor to ensure that the processing, as actually performed using equipment, materials, and personnel in the processing facility, achieves the desired result. This requires verification and/or validation and routine monitoring of the process.

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