




A New Standard of Care for OSA Treatment

Inspire V™ System



Designed for a simplified procedure
& enhanced patient comfort.

Inspire therapy has transformed the lives of over **100,000 patients worldwide**¹.

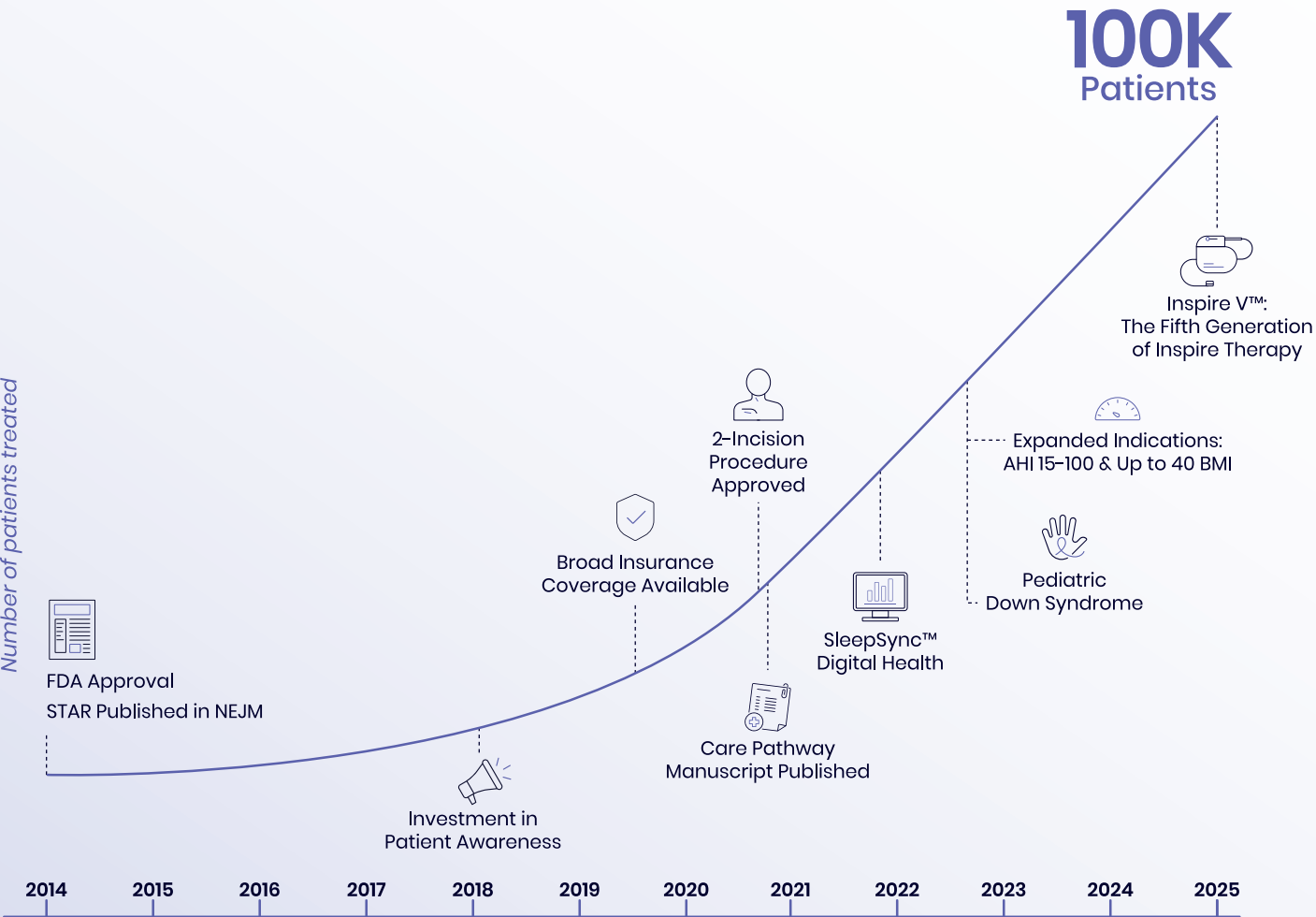

¹ Data on file.

OSA affects over
30 MILLION AMERICANS

When left untreated OSA can lead to serious health issues.

INCREASED RISK OF:

Stroke • Sudden Cardiac Death • Risk of Cardiovascular Mortality²



Proven Outcomes with Inspire Therapy

Inspire therapy is a standard of care for people who struggle with OSA and aren't able to benefit from CPAP.



² Jean-Louis G et al. Obstructive sleep apnea and cardiovascular disease: role of the metabolic syndrome and its components. J Clin Sleep Med 2008;4(3):261-272.

³ Rotenberg et al. Trends in CPAP adherence over twenty years of data collection: a flattened curve. Journal of Otolaryngology Head and Neck Surgery 2016; 45:43

⁴ Data on file.

⁵ Woodson et al. Upper airway stimulation for obstructive sleep apnea: 5-year outcomes. OTO-HNS 2018

Simplified Procedure

One Less Step

With a **new integrated respiratory sensor**, Inspire V eliminates the need to place a separate sensing lead.

Designed for patient comfort

Advanced Customization

Tailored therapy options enable more gradual changes in stimulation levels allowing more patients to achieve the full benefits of OSA therapy.




NEW integrated sensing

Proprietary technology allows for respiratory sensing to occur inside of the Inspire V implant.

NEW stimulation step sizes


provides precise therapy stimulation options.

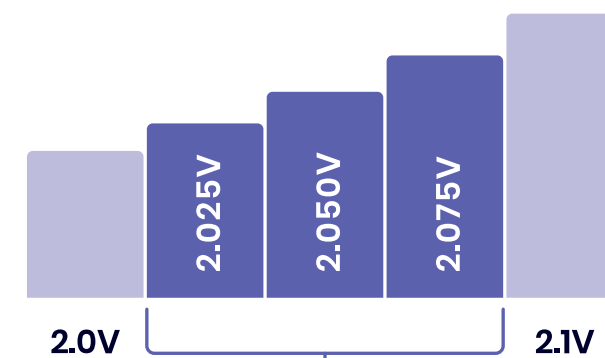

Program patient therapy levels at half or quarter step sizes with Inspire V when needed.

NEW ramp option

designed to gently initiate therapy, without disturbing sleep.

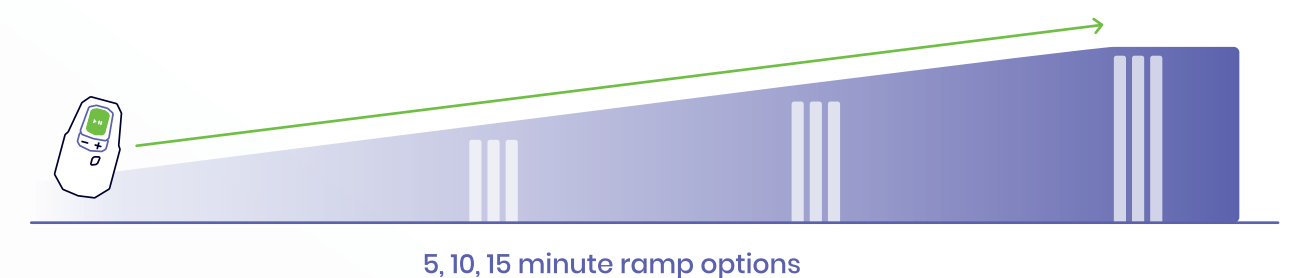

	Inspire V	Inspire IV
Simplified Procedure		
Integrated Respiratory Sensor	✓	
Designed for Patient Comfort		
Ramp Option for Therapy Onset	✓	
Advanced Options for Stimulation Step Sizes	✓	
Consistent Proven Benefits		
Patient Controlled Therapy	✓	✓
Synchronized Stimulation Algorithm	✓	✓
Inspire SleepSync™ Digital Health Platform	✓	✓
Inspire® App	✓	✓

Inspire therapy is indicated for patients who are*:

- ✓ Diagnosed with moderate to severe OSA (AHI 15-100)
- ✓ Unable to use or get consistent benefit from CPAP
- ✓ Not significantly overweight (BMI ≤ 40)
- ✓ Over 18 years of age



MRI Conditional Full body 1.5T MRI scans can be performed following certain guidelines**

New System

Model Number	Description
3150	Implantable Pulse Generator
4063	Stimulation Lead
2580	Sleep Remote



*Visit www.inspiresleep.com for more information

**For full MRI requirements, warnings, and precautions, please view the "MRI Guidelines for Inspire Therapy" manual at manuals.inspiresleep.com

Inspire is not for everyone. It is a surgically implanted system that is intended to treat obstructive sleep apnea in patients who are not effectively treated by, or able to tolerate CPAP. Talk to your patients about risks, benefits and expectations associated with Inspire. Risks associated with the surgical implant procedure may include infection and temporary tongue weakness. In rare cases tongue paresis and atrophy may occur. Some patients may require post implant adjustments to the system's settings in order to improve effectiveness and ease any initial discomfort they may experience. Important safety information and product manuals can be found at inspiresleep.com/safety-information/ or call 1-844-OSA-HELP

Inspire and the cloud design are registered trademarks of Inspire Medical Systems, Inc. © 2025 Inspire Medical Systems, Inc. All rights reserved. 802-222-001, Rev E.