

9931 - Halo Sectional Matrix System

Original Bands with Instruments Kit

- 1 x Each *Universal Ring, Short Ring*
- 25 x Each *Original Matrices*
- 3.5 mm, 4.5 mm, 5.5 mm, 6.5 mm, and 7.5 mm
- 25 x Each *Wedges Extra Small, Small, Medium, and Large*
- 1 x *Halo Carousel* 1 x *Tweezers* 1 x *Forceps*

9932 - Halo Sectional Matrix System

Original Bands Kit

- 1 x Each *Universal Ring, Short Ring*
- 25 x Each *Original Matrices*
- 3.5 mm, 4.5 mm, 5.5 mm, 6.5 mm, and 7.5 mm
- 25 x Each *Wedges Extra Small, Small, Medium, and Large*
- 1 x *Halo Carousel*

9933 - Halo Sectional Matrix System

Original Bands Mini Kit

- 1 x Each *Universal Ring, Short Ring*
- 5 x Each *Original Matrices*
- 3.5 mm, 4.5 mm, 5.5 mm, 6.5 mm, and 7.5 mm
- 5 x Each *Wedges Extra Small, Small, Medium, and Large*
- 1 x *Halo Carousel*

9935 - Halo Sectional Matrix System Firm Nonstick Bands Kit

- 1 x Each *Universal Ring, Short Ring*
- 20 x Each *Firm Nonstick Matrices*
- 3.5 mm, 4.5 mm, 5.5 mm, 6.5 mm, and 7.5 mm
- 25 x Each *Wedges Extra Small, Small, Medium, and Large*
- 1 x *Halo Carousel*

9936 - Halo Sectional Matrix System Firm Bands Kit

- 1 x Each *Universal Ring, Short Ring*
- 25 x Each *Firm Matrices*
- 3.5 mm, 4.5 mm, 5.5 mm, 6.5 mm, and 7.5 mm
- 25 x Each *Wedges Extra Small, Small, Medium, and Large*
- 1 x *Halo Carousel*



Halo Original Matrices	50pk	100pk
3.5 mm	5448	5449
4.5 mm	5450	5451
5.5 mm	5452	5453
6.5 mm	5454	5455
7.5 mm	5456	5457



Halo Firm Matrices	50pk	100pk
3.5 mm	5059	—
4.5 mm	5062	5063
5.5 mm	5064	5065
6.5 mm	5066	5067
7.5 mm	5068	—



Halo Firm Nonstick Matrices	50pk	100pk
3.5 mm	5049	—
4.5 mm	5051	5052
5.5 mm	5053	5054
6.5 mm	5055	5056
7.5 mm	5057	—

Halo Wedges	100pk
Extra Small	5884
Small	5042
Medium	5043
Large	5044

NEW

Size	1pk	2pk
Universal	5008	5009
Short	5882	5883

NEW

5047 - Halo Forceps 1pk



5048 - Halo Tweezers 1pk



NEW



SHORT RING

NEW



EXTRA SMALL WEDGE

Beautiful Results in Less Time

The easy-to-use Halo sectional matrix system allows you to create beautiful, anatomically contoured composite restorations in less time.¹



HALO™

SECTIONAL MATRIX SYSTEM

The Halo system produces anatomically shaped contours so you can spend less time shaping and finishing your restoration.¹ And it provides consistently tight contacts, helping to ensure that your restoration will have the longevity both you and your patient expect.

Maximum tooth separation is achieved through the force of the nitinol ring. This ideal separation allows you to restore a single tooth or back-to-back restorations with ease.

The anatomically contoured design of the Halo ring, along with the adaptation provided by the wedges, allows the band to fully adapt to the prep from the gingival margin to the marginal ridge and helps it maintain its shape even in large preparations, helping you to consistently achieve anatomically shaped restorations.

When your restorations have ideal contacts and anatomical curvature, you will see a reduction in the time spent shaping and finishing.¹ The Halo system is designed to do that with every restoration.



The Halo sectional matrix system includes nitinol rings, stainless steel matrix bands, and wedges. Each of the components were designed to work together to create consistent, anatomically shaped restorations.

NEW



NITINOL RINGS

- New Short ring provides the same separation with a better fit on short clinical crowns
- Super-elastic nitinol metal maintains force during procedures and reduces cyclic fatigue
- Rings will last over 1,000 uses¹
- Glass-filled nylon tines won't easily break or collapse into the prep and create under contoured restorations
- Ring contours secure the band in a natural and anatomical shape, helping to eliminate flash and reducing the amount of finishing needed
- Stackable design allows for use with MODs and other Class II restorations
- Provides ideal separation for back-to-back restorations



MATRIX BANDS

- Anatomically shaped for ideal contacts
- Curve at marginal ridge creates ideal occlusal embrasure, reducing finishing
- Tweezer holes for easy placement
- Bendable tab allows for easy placement and removal with increased visibility
- Optional coating for easy identification
- Original bands allow for more adaptability and burnishing, while Firm Bands resist deformation and are ideal for tight interproximal spaces

NEW



WEDGES

- New Extra Small wedge
- Firm wedge is less traumatic to the papillae
- Hollow design makes wedges easy to place and allows wedges to be stacked when multiple wedges are needed
- Easy to distinguish colors help identify sizes of wedges
- Collapsible center for anatomical adaptation of the band



CAROUSEL DISPENSER

- Stay organized with unique, easy-to-use carousel dispenser
- No more spilling and scattering of bands and wedges after an accidental drop, saving time and frustration
- Can be easily disinfected with a surface disinfectant

Case photos courtesy of Dr. Jaleena Fischer-Jessop, Chief Clinical Officer of Ultradent Products, Inc.




Pre-op



Preparation with Halo matrix system



Post-op with Transcend composite Universal Body shade

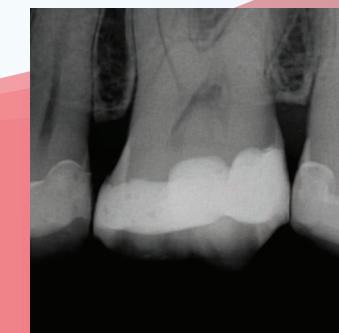

Post-op radiograph

Case photos courtesy of Dr. James Wright, Clinical Advisor of Ultradent Products, Inc



Pre-op



Preparation with Halo system



Post-op

Case photos courtesy of Dr. James Wright, Clinical Advisor of Ultradent Products, Inc



Pre-op



Preparation with Halo matrix system



Post-op